

Алматы (7273)495-231  
Ангарск (3955)60-70-56  
Архангельск (8182)63-90-72  
Астрахань (8512)99-46-04  
Барнаул (3852)73-04-60  
Белгород (4722)40-23-64  
Благовещенск (4162)22-76-07  
Брянск (4832)59-03-52  
Владивосток (423)249-28-31  
Владикавказ (8672)28-90-48  
Владимир (4922) 49-43-18  
Волгоград (844)278-03-48  
Вологда (8172)26-41-59  
Воронеж (473)204-51-73  
Екатеринбург (343)384-55-89

Ижевск (3412)26-03-58  
Иваново (4932)77-34-06  
Иркутск (395)279-98-46  
Казань (843)206-01-48  
Калининград (4012)72-03-81  
Калуга (4842)92-23-67  
Кемерово (3842)65-04-62  
Киров (8332)68-02-04  
Коломна (4966)23-41-49  
Кострома (4942)77-07-48  
Краснодар (861)203-40-90  
Красноярск (391)204-63-61  
Курск (4712)77-13-04  
Курган (3522)50-90-47  
Липецк (4742)52-20-81

Киргизия (996)312-96-26-47

Магнитогорск (3519)55-03-13  
Москва (495)268-04-70  
Мурманск (8152)59-64-93  
Набережные Челны (8552)20-53-41  
Нижний Новгород (831)429-08-12  
Новокузнецк (3843)20-46-81  
Ноябрьск (3496)41-32-12  
Новосибирск (383)227-86-73  
Омск (3812)21-46-40  
Орел (4862)44-53-42  
Оренбург (3532)37-68-04  
Пенза (8412)22-31-16  
Петрозаводск (8142)55-98-37  
Псков (8112)59-10-37  
Пермь (342)205-81-47

Россия (495)268-04-70


Ростов-на-Дону (863)308-18-15  
Рязань (4912)46-61-64  
Самара (846)206-03-16  
Саранск (8342)22-96-24  
Санкт-Петербург (812)309-46-40  
Саратов (845)249-38-78  
Севастополь (8692)22-31-93  
Симферополь (3652)67-13-56  
Смоленск (4812)29-41-54  
Сочи (862)225-72-31  
Ставрополь (8652)20-65-13  
Сыктывкар (8212)25-95-17  
Сургут (3462)77-98-35  
Тамбов (4752)50-40-97  
Тверь (4822)63-31-35

Казахстан (772)734-952-31

Тольяти (8482)63-91-07  
Томск (3822)98-41-53  
Тула (4872)33-79-87  
Тюмень (3452)66-21-18  
Улан-Удэ (3012)59-97-51  
Ульяновск (8422)24-23-59  
Уфа (347)229-48-12  
Хабаровск (4212)92-98-04  
Чебоксары (8352)28-53-07  
Челябинск (351)202-03-61  
Череповец (8202)49-02-64  
Чита (3022)38-34-83  
Якутск (4112)23-90-97  
Ярославль (4852)69-52-93

<https://integra.nt-rt.ru> || [iag@nt-rt.ru](mailto:iag@nt-rt.ru)

# Программное обеспечение для автоматизации пипеток VIALAB

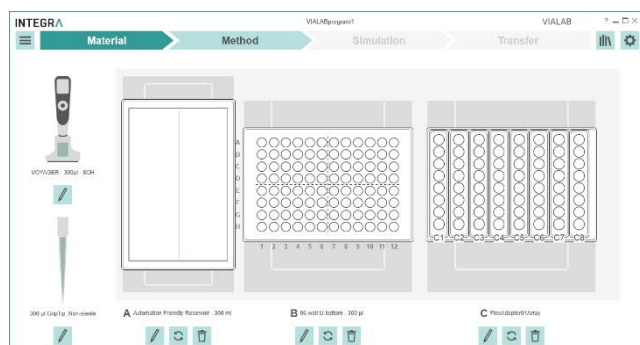

## VIALAB Operating Instructions

### Document Change History

<b>Vers.</b>	<b>Status:</b>	<b>Name:</b>	<b>Date:</b>	<b>Change:</b>
V00	Draft	ULe, SSt	07-12-2017	Document issued
V01	Released	SSt	15-03-2018	Document released
V02	Released	SSt	10-10-2018	download information added, troubleshooting "Program Transfer To Pipette" added
V03	Released	ULe, SSt	09-01-2019	New Software V 2.0
V04	Released	SSt	13-02-2019	System Requirements for VIALAB 2.0 changed
V05	Released	SSt	04-03-2020	New Software V 2.1 – new features
V06	Released	MFr	21-04-2020	Update Pipette Communication Module
V07	Released	JBa2	13-08-2020	Typos removed
V08	Released	ULe, JBa2, SSt	08-04-2022	New Software V 3.0
V09	Released	SSt	05-12-2022	System Requirements updated

## Table of Contents

<b>1</b>	<b>Introduction</b>	<b>4</b>
1.1	Purpose	4
1.2	Scope	4
1.3	Third party licenses	4
<b>2</b>	<b>System requirements</b>	<b>5</b>
2.1	Personal computer / laptop	5
2.2	Compatibility	5
2.3	Useful accessories	5
<b>3</b>	<b>Installation of VIALAB</b>	<b>6</b>
3.1	Preparation	6
3.2	First time installation	6
3.3	Running VIALAB without installation	6
<b>4</b>	<b>Settings</b>	<b>7</b>
4.1	VIALAB settings and information	7
4.2	Update VIALAB	7
<b>5</b>	<b>Main menu</b>	<b>8</b>
5.1	Information and saving options	8
5.1.1	Program Settings	9
<b>6</b>	<b>Programming</b>	<b>10</b>
6.1	Overview	10
6.2	Material	12
6.3	Method	14
6.4	Method steps only for D-ONE pipetting module	23
6.4.1	Importing worklists	23
6.4.2	Hit Picking	27
6.4.3	Normalization	29
6.4.4	Normalization Worklist	31
6.4.5	Master Mix	33
6.4.6	Worklist	34
6.5	Simulation	35
6.6	Transfer	36
6.6.1	Connect pipette or pipetting module to the PC	36
6.6.2	Transfer	36
6.7	Run a program directly from VIALAB	37
<b>7</b>	<b>Labware library</b>	<b>38</b>
7.1	Reservoirs	38
7.2	Plates	40
7.3	Tube Rack	42
7.4	COMBI Base	43
7.5	COMBI System	44
7.5.1	Slanted Plate Holder	44
<b>8</b>	<b>Trouble Shooting</b>	<b>46</b>

## 2 System requirements

### 2.1 Personal computer / laptop

CPU: INTEL or AMD (not older than 5 years)  
USB port: 2.0 / 3.0  
Bluetooth: optional  
Memory: > 4 GB  
Internet access: recommended during installation, required for updates  
Operating system: MS Windows 8.1, MS Windows 10.x, MS Windows 11.x,  
Framework: Microsoft .NET 4.7.2 (or later)  
Resolution: 1920 x 1080, Scaling/Zoom factor 100 %

### 2.2 Compatibility

Pipette:	Serial No.:	Pipette Firmware connectable to			
		VIALAB V 1.0.x	VIALAB V 2.0.x	VIALAB V 2.1.x	VIALAB V 3.0.x
VIAFLO electronic multichannel pipette	7000000 and higher	4.12 and higher	4.20 and higher	4.21 and higher	5.04 and higher
VOYAGER electronic adjustable tip spacing pipette	7000000 and higher	4.12 and higher	4.20 and higher	4.21 and higher	5.04 and higher
D-ONE single channel pipetting module	7000000 and higher	-	-	-	5.04 and higher

### 2.3 Useful accessories


Part No.	Name	Description
4211	Charging/communication stand for 1 pipette, incl. mains adapter	Connect the stand to the USB port of your PC and place the pipette/D-ONE on the stand. The stand can be used to charge the battery and enables the communication between the PC and the pipette/D-ONE.
4221	Pipette communication module	For wireless communication between pipette/D-ONE and PC or ASSIST PLUS.
4549	Charging cable for pipette on ASSIST PLUS	For charging of pipette/D-ONE on ASSIST PLUS.

## 4 Settings

### 4.1 VIALAB settings and information



Click on the Settings button in the top right corner. Use the VIALAB Settings to define your preferences and get additional information.

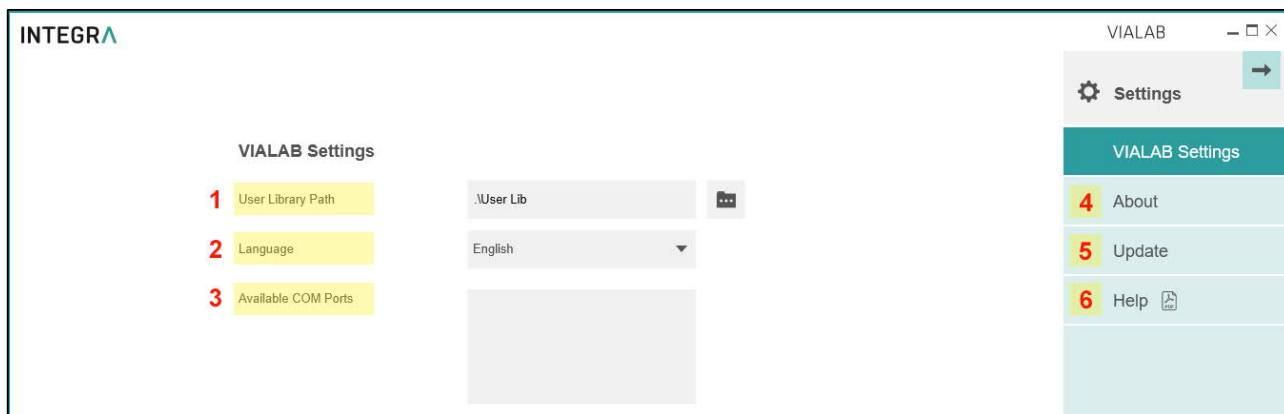

**1 User Library Path:** define the storage path of your labware library.

**2 Language:** choose your preferred language.

**3 Available COM Ports:** shows all available communication ports.


**4 About:** information about the VIALAB software version and the INTEGRA contact data

**5 Update:** check for new updates, see 4.2.

**6 Help:** download the latest VIALAB operating instructions.

### 4.2 Update VIALAB

If an update for VIALAB is available and you are connected to the internet you can check for new updates.

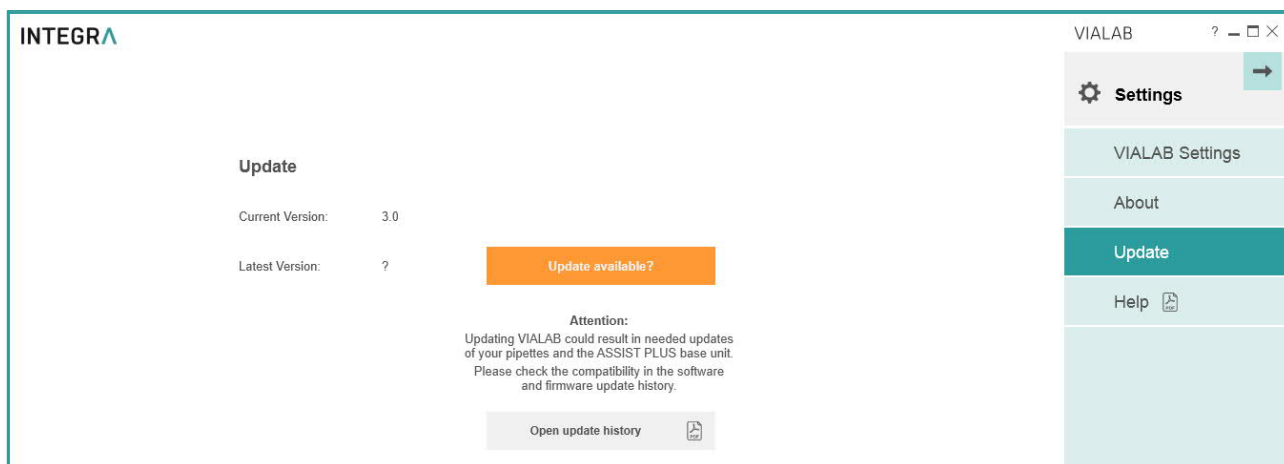

Click on Open update history to download information about software versions of VIALAB and VIALINK as well as firmware information of VIAFLO/VOYAGER pipettes, D-ONE module and ASSIST / ASSIST PLUS.

## 5 Main menu

### 5.1 Information and saving options



Within a program, click on the Menu button in the top left corner to open the main menu. Here information about the program is displayed and saving options are listed.

The screenshot shows the VIALAB SerialDilution program interface. The top bar includes the INTEGRAL logo, the program name 'SerialDilution', and the VIALAB logo. A sidebar on the left contains a menu with 10 numbered options: 1. Information, 2. Program Settings, 3. Save, 4. Save As, 5. Report, 6. Close, 7. New, 8. Open, 9. Sample Programs, and 10. Exit. The main area displays program details for 'SerialDilution', including a description, material, and pipette information. Below this is a 'Program overview' section showing a deck layout with pipette, tip, and deck content details.

Pipette	Tip	Deck Content
VOYAGER - 300ul - 8CH	300 ul GripTip, Sterile, Filter, Low retention	Initial Volumes, Transfer

**1 Information:** change the Program name and add a description, if desired.

**2 Program Settings:** define the clearance height and offsets of each deck position for your program.

**3 Save:** save your program.

**4 Save As:** save you program under a new name.

**5 Report:** generate a printable PDF-report of your program, e. g. for standard operating procedures.

**6 Close:** close your program.

**7 New:** create a new program.

**8 Open:** open an existing program, recently used programs are listed.

**9 Sample Programs:** select and open a sample program.


**10 Exit:** exit VIALAB.



**Note:** You should save your current VIALAB program on your computer to have a backup or if you need to transfer the program to more than one pipette.

### 5.1.1 Program Settings


The following settings apply to each VIALAB program individually.

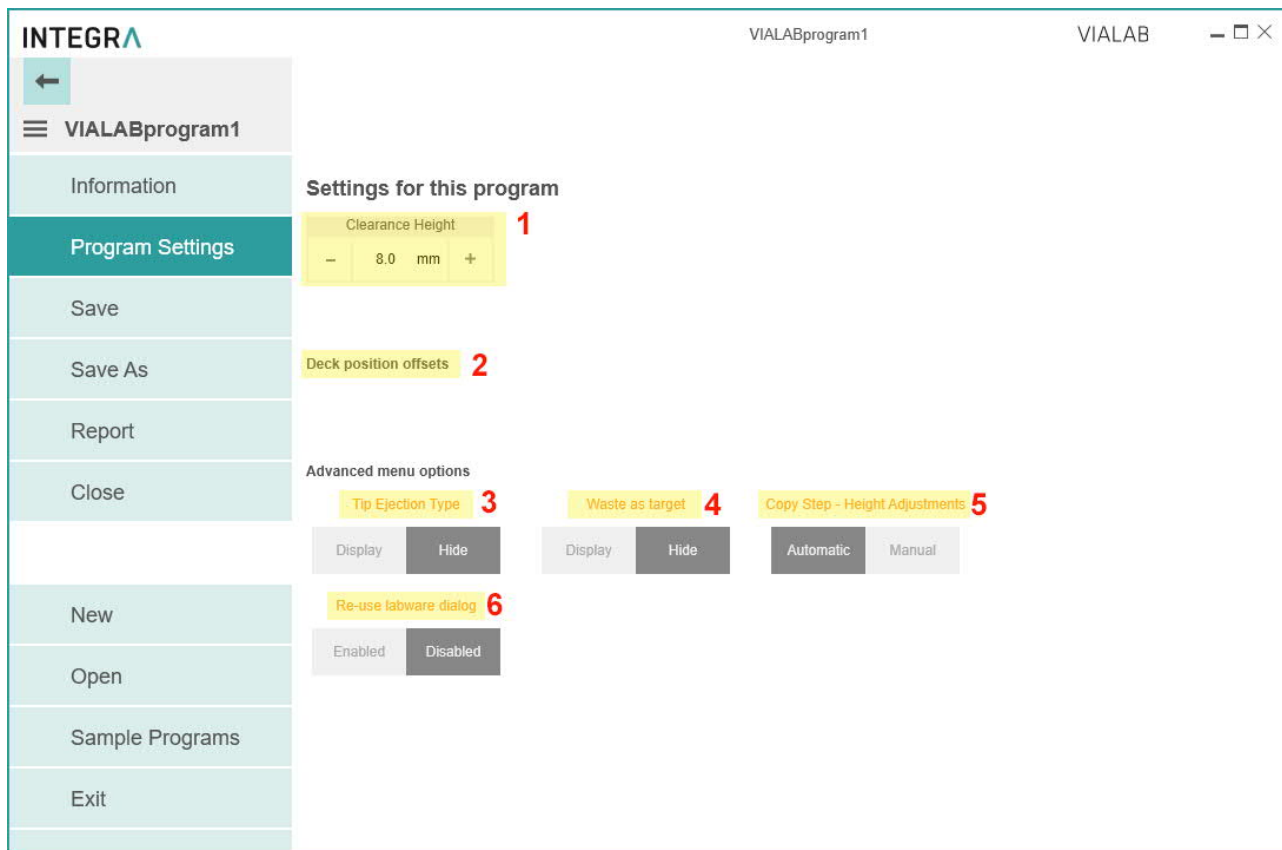

- 1 Clearance Height:** distance over the highest point of the labware at which the pipette moves over the labware.
- 2 Deck Position Offsets:** deck positions A, B or C can be shifted in X, Y and Z axis for max. 2 mm.
- 3 Tip Ejection Type:** choose whether to display or hide an additional option for the tip ejection. If activated in the tab “Tip Change” the tips can be ejected either at the ejector ledge (**Ledge**) or directly in the waste bin without touching the ejector ledge (**Waste**) during tip ejection.
- 4 Waste as target:** choose whether to display or hide the option to dispense directly into the waste bin. If activated the waste bin can be selected as target position in a Multi Aspirate step or in Serial Dilution step as a waste position for the last aspiration.
- 5 Copy Step – Height Adjustments:** when copying steps, the liquid level changes for the copied step. If tip travel is activated in a step, the heights are automatically adjusted relative to the new liquid level (Automatic). Alternatively, the heights of a copied step aren’t changed and identical to the original step (Manual).
- 6 Re-use Labware Dialog:** it is helpful to enable this function when working with more than one labware layout (see section 6.2). If you (accidentally) delete a labware in one of the following layouts, clicking on its position opens a dialog from which the recently used labware can be re-used.

## 6 Programming

### 6.1 Overview

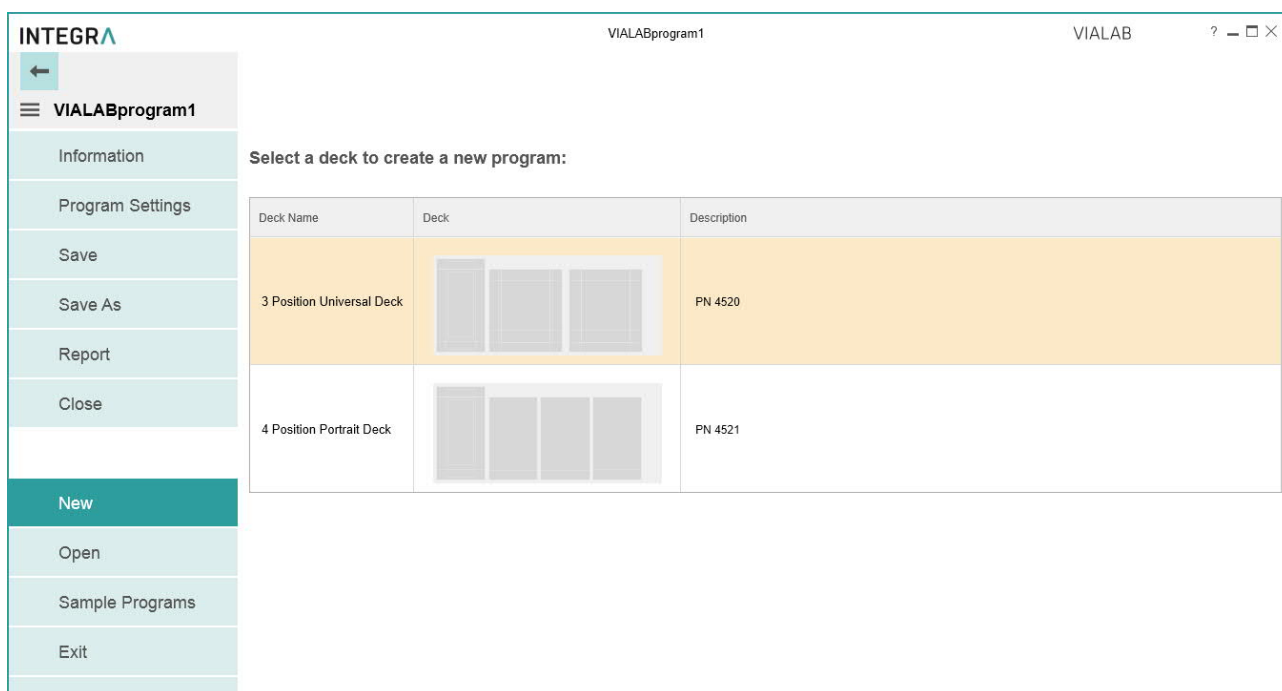

After starting VIALAB, you can create a new pipetting program, open an existing or recently used program or get started with a sample program.


To create a new program, select first the deck you want to use, e.g. the 3 Position Universal Deck with one reservoir position and two spring loaded plate positions.



**Note:** The D-ONE pipetting module requires its own tip deck (PN 4535) which need to be installed on the ASSIST PLUS.



The program menu will be displayed:

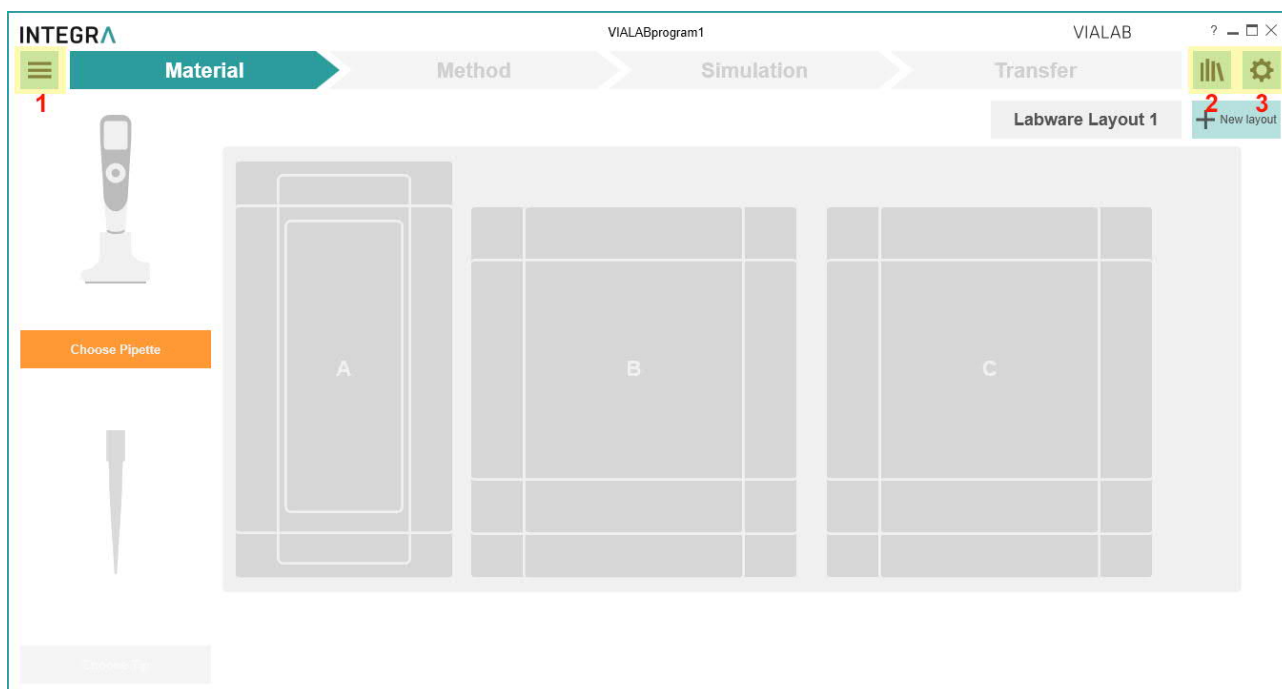


- 1 Main Menu:** information and saving options of the defined program (see chapter 5).
- 2 Library:** labware library listing available decks, pipettes, tips, reservoirs, plates, tube racks, and COMBI systems with COMBI bases which can be used for a program. You can also define new reservoirs, plates, tube racks and COMBI systems. Alternatively, import already defined labware files to the labware library.
- 3 Gear Wheel:** settings and updates of VIALAB.

The table below shows the different sections for creating a program. After one section is completed, the next can be defined.


Section	Description
Material	Choose the pipette with tips and the labware on the deck, e.g. the reservoir used as source container and the kind of well plates with orientation.
Method	Create your individual pipetting program
Simulation	Review a simulation of all steps in your pipetting program.
Transfer	Transfer your VIALAB program to your pipette and run it from your pipette. Alternatively connect the ASSIST PLUS via USB to the PC, connect the Pipette and start a program from VIALAB directly.

## 6.2 Material

First click on the orange square to choose any multichannel VIAFLO/VOYAGER pipette or D-ONE module and afterwards select the GripTips. This enables the selection of the other labware used in the program.

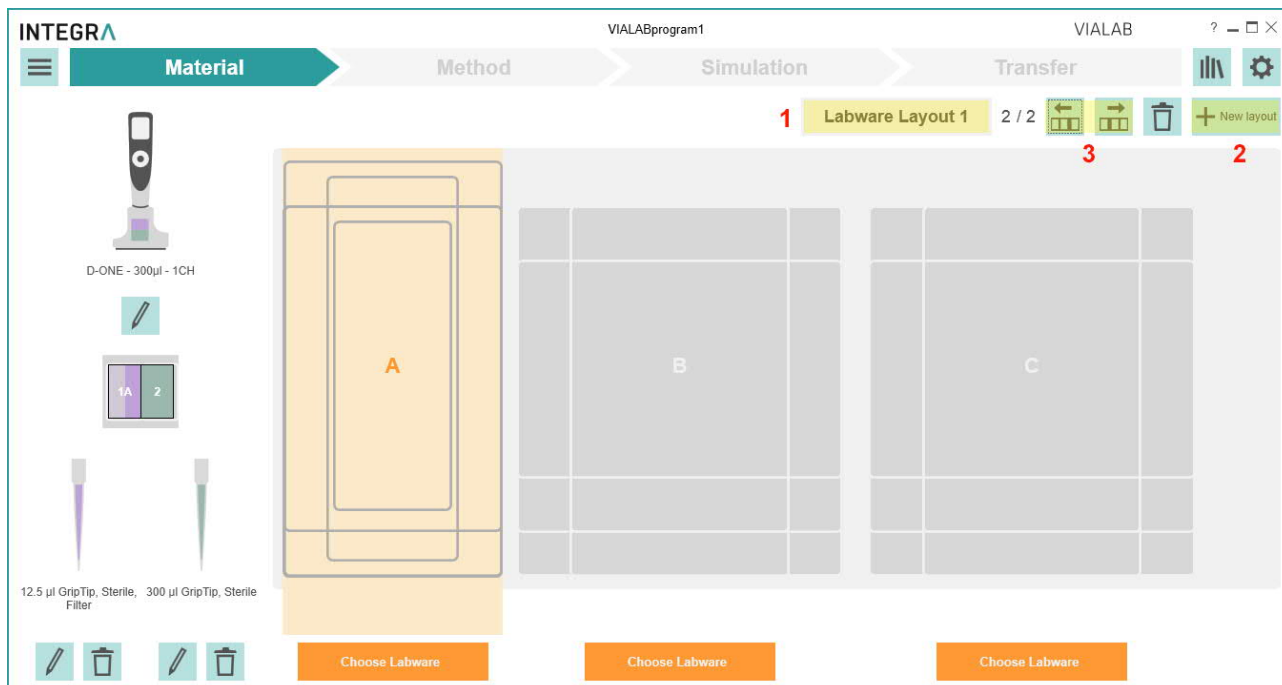


Click on an orange button or directly on the deck position to choose the predefined labware on all required positions. Click on the Labware Layout button (1) to name the deck. To create an additional labware layout, click on New layout (2). Use the arrow keys to switch between all available labware layouts (3).

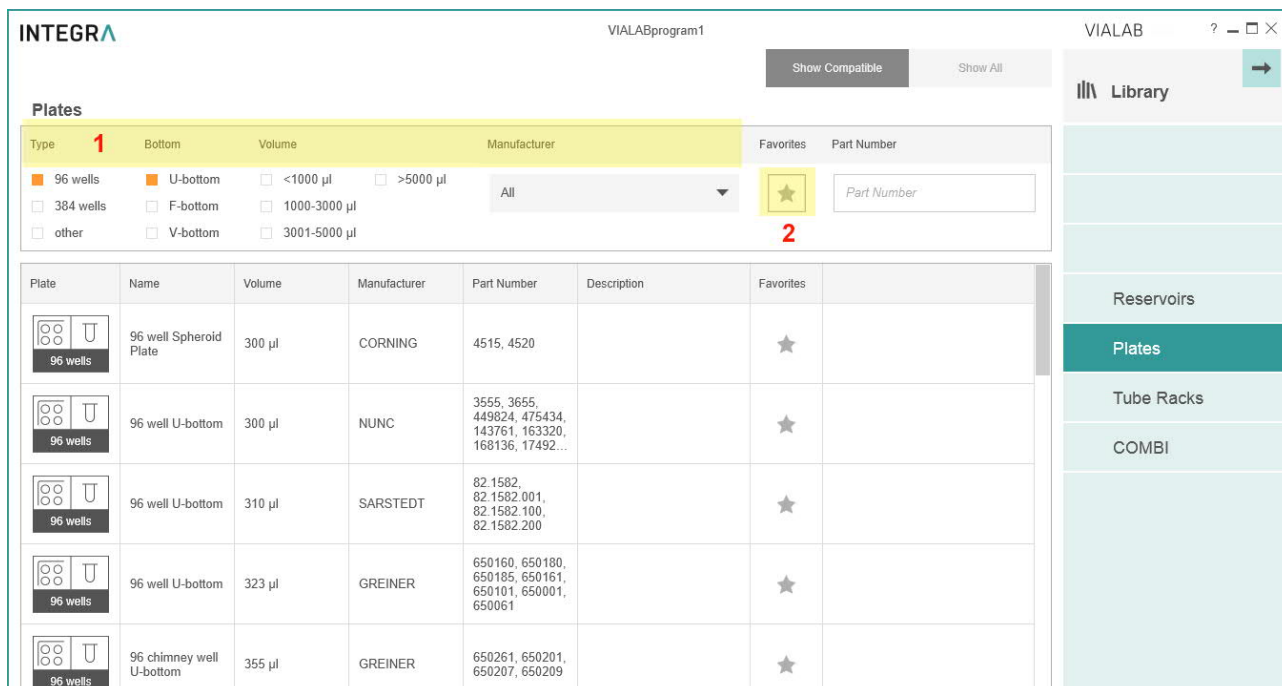

Select the desired properties or manufacturer (1) to quickly show reservoirs/plates/tube racks/COMBI Systems which are compatible with the selected pipette and deck position. Click on the star (2) to show favorites only (to be defined in the labware library).

Scroll to one reservoir, plate, tube rack or COMBI System and select it by clicking on it.

Labware	Description
Reservoirs	Select an INTEGRA reservoir or define a new reservoir, see Library icon.
Plates	Select one of predefined, commercially available plates. Search for any plate or enter the part number directly, if known. Alternatively, define the dimensions of a new plate, see Library icon.
Tube Racks	Select a tube rack or define a new one, see Library icon.
COMBI System	To add a COMBI System, select a COMBI BASE and fill each selected column with modules, e.g. plates, tube racks or reservoirs and click on Apply COMBI System, see also section 7.5.

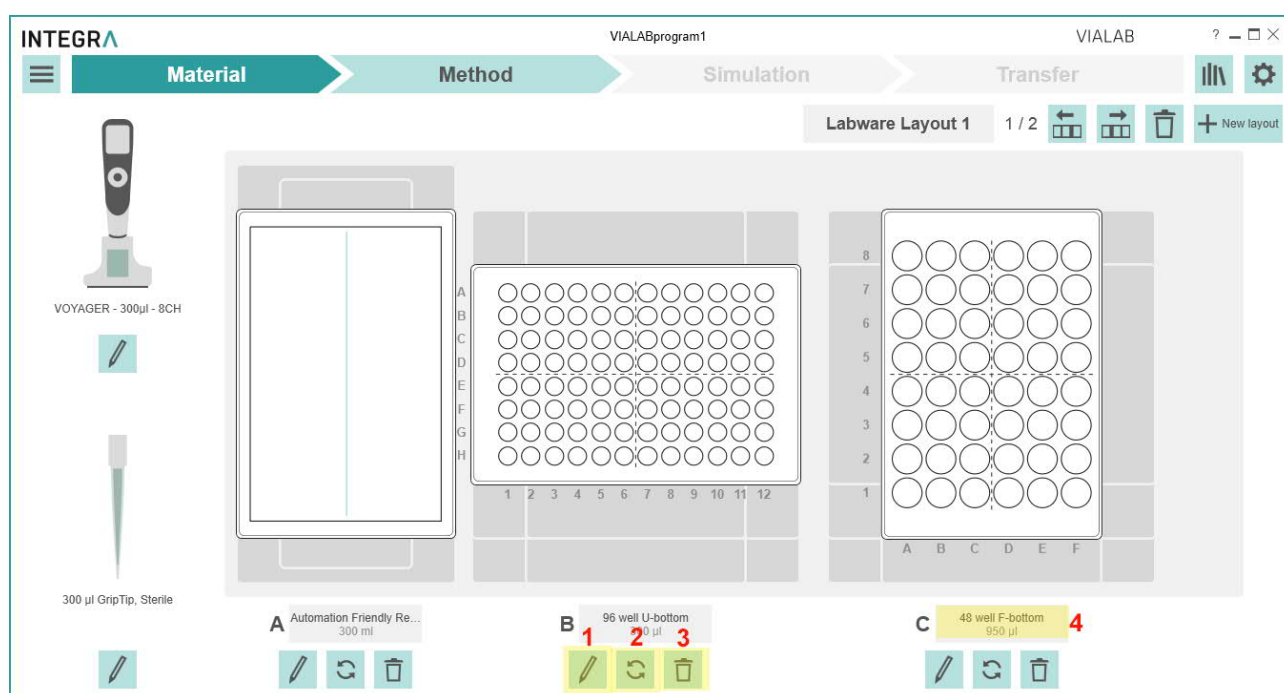

**1 Edit:** change the chosen labware.

**2 Change orientation:** rotate counterclockwise by 90° (all 4 possible orientations).

**3 Delete:** delete the selected labware on this position.


**4 Rename:** name the chosen labware.

**Note:** changing the material of an already defined program, e.g. pipette or plate orientation, will delete all subsequent program steps in the “Method” section, if they are not compatible with the previous settings.

## 6.3 Method

Your pipetting program is composed of several steps.

When you open the section Method the first time, the first step “Initial Volumes” is already preset, and you are prompted to set the volumes of the labware which are already filled.

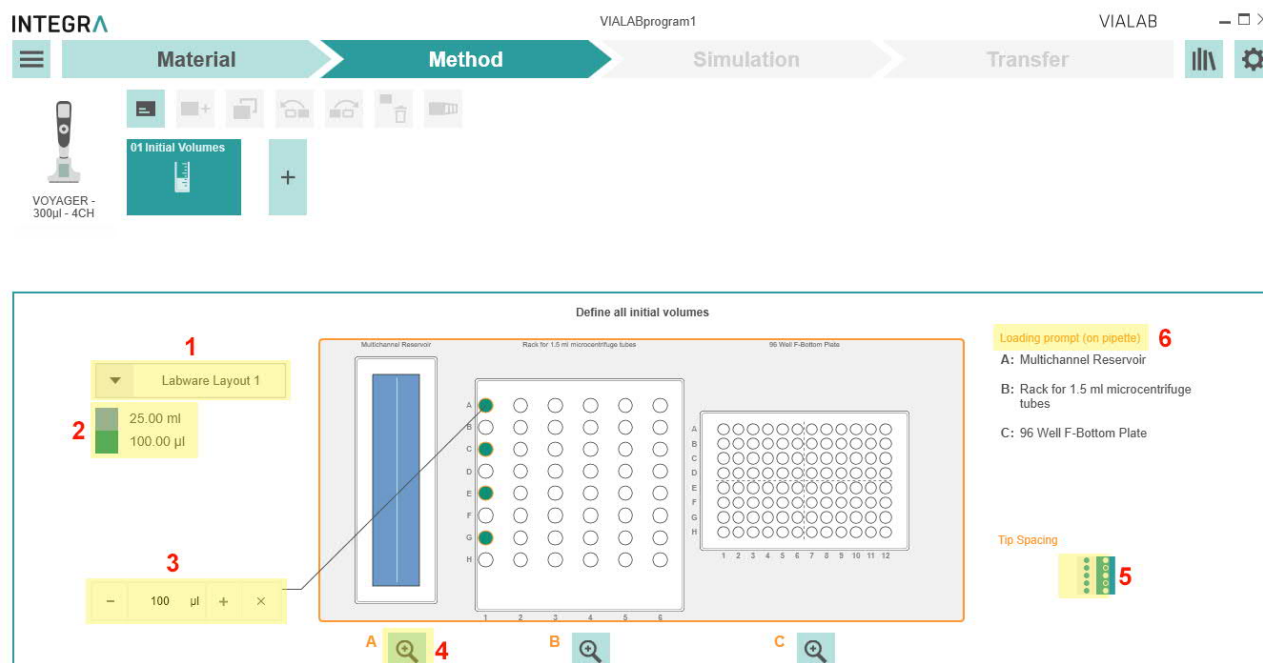

- 1 Choose labware layout:** if more than one is available, choose the corresponding layout.
- 2 Color legend:** shows the pre-fill (initial) volume of the labware on position A, B or C. Each color represents one pre-fill volume.
- 3 Pre-fill volume:** click on the vessel which is already filled at the start of your program and enter the pre-fill volume.
- 4 Zoom:** displays a close-up view of one deck position.
- 5 Tip spacing:** choose a different tip spacing to change the filling pattern of your plate e.g. as shown with plate on position C (only for VOYAGER tip spacing pipettes).
- 6 Loading prompt (on pipette):** shows the used labware inclusive deck position of every available labware layout within the current VIALAB program.


Select the wells which are prefilled at the start of your program via click or drag and drop. In the pop-up window you type in the dedicated volume in the given unit to assign a volume.

To view or edit the assigned volume select wells and the current volume is displayed. If needed, edit the volume, or delete it.







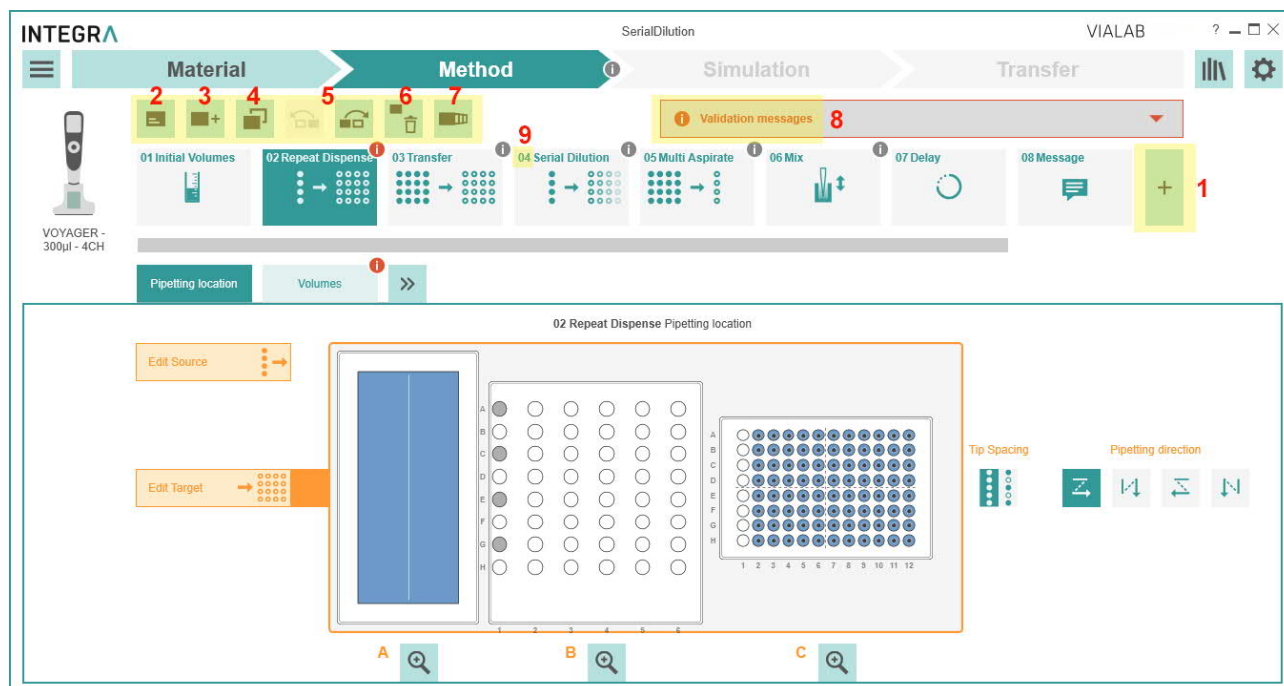





**Note:** If you're using D-One pipetting modules, this step can be skipped as D-One features liquid level detection.

To add further steps, click on a plus square (1, 3). To move, copy or delete steps click on the corresponding icon (4, 5, 6). The steps are processed sequentially according to the step number (9).





- 1 Plus square:** add a new step to the end of the method.
- 2 Add Description:** add a description for the current step. The first 15 characters are displayed at the bottom of the step icon.
- 3 Add Step:** add a new step after a selected step to the method.
- 4 Copy Step:** copy selected step to the end.
- 5 Move left/right:** move the selected step forwards or backwards.
- 6 Delete Step:** delete selected step.
- 7 Multi-Selection:** click on this button to select multiple steps in parallel to copy or delete them.
- 8 Validation messages:** click on this field for programming error details.
- 9 Step number:** the step number indicates the step sequence in the program.

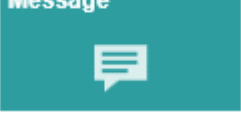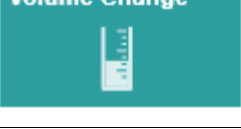Following steps are available for VIAFLO/VOYAGER pipettes and D-ONE modules:


Symbol	Step	Description
<b>Transfer</b> 	Transfer	Transfers an equal number of samples from a source to target (N:N), e.g. plate copy.
<b>Repeat Dispense</b> 	Repeat Dispense	Dispenses from one source into multiple targets (1:N), e.g. reagent addition
<b>Serial Dilution</b> 	Serial Dilution	Aspirates sample from a source and performs a serial dilution in the defined number of wells.
<b>Multi Aspirate</b> 	Multi Aspirate	Collects several samples from one source and transfers into one target (N:1).
<b>Mix</b> 	Mix	Mixes a sample by pipetting up and down in the selected vessel.
<b>Delay</b> 	Delay	Defines a pause between the last and the next step in seconds.
<b>Message</b> 	Message	Pauses the program and displays a message. Three lines with 12 characters each are available.
<b>Volume Change</b> 	Volume Change	Enter the filling volume of the vessels e.g. after manual filling steps within a method. The volume defined in this step is an absolute volume. It does not consider previous liquid handling steps.
<b>AUX Port</b> 	AUX Port	Turning devices on and off which are connected to the auxiliary port of the ASSIST PLUS.
<b>Labware Change</b> 	Labware Change	Change labware items on the deck within the same program. Switch between the labware layouts defined under Materials.

For D-ONE pipetting module additional steps are available, see 6.4.


Within one step, adjust the following parameters:

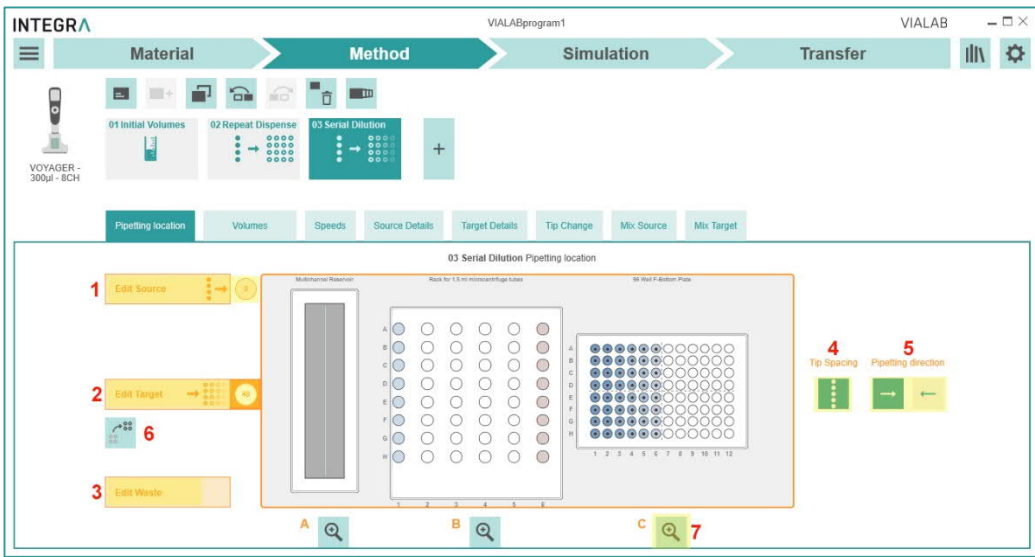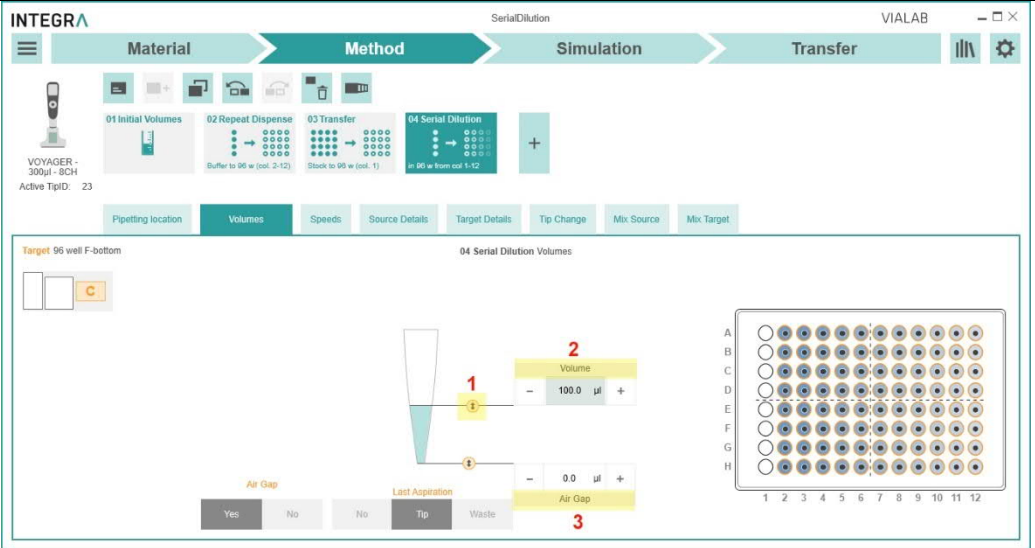
Parameter	Description
Pipetting location	<ul style="list-style-type: none"> <li><b>Edit Source (1):</b> click on the button and then select the wells (for plates on upmost well of column) where the liquid should be aspirated from. It will be colored. The number of source wells is displayed next to it, e.g. 8.</li> <li><b>Edit Target (2):</b> click on the button and then select the wells (for plates on upmost well of column) where the liquid should be dispensed in. It will be marked with a black dot in the middle of the well. The number of target wells is displayed next to it, e.g. 48.</li> <li><b>Edit Waste (3):</b> only for step “Serial Dilution” and if under Volumes “Last aspiration back to “Waste” was selected. Click on the button to select waste vessel location.</li> </ul>  <p>If needed, adapt the tip spacing (4) or the pipetting direction (5). For changing the selected wells deselect the chosen ones or move all selected wells together within one plate (6). For a detailed view of the plates, use the zoom-in function (7).</p>
Volumes	 <p>Move the slider up or down (1), type in the volume (2) or use the “+” and “-” buttons to set the desired <b>Volume(s)</b> for the current step. An <b>Air Gap</b> (3) can be set to aspirate a defined air volume after an aspiration of liquid.</p>

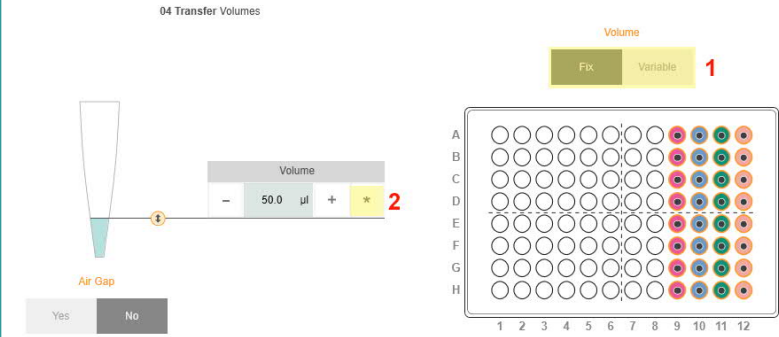

Volumes	Step Type	Volume description
	Transfer	<p><b>Volume:</b> define whether the volume is variable or fix for all transfers (<b>1</b>). If variable, select each individual column and assign a volume. If the calculated volume exceeds the tip volume, the volume is dispensed multiple times with the same tip, indicated by a multiplier (*,<b>2</b>).</p> 
	Repeat Dispense	<p><b>Post-Dispense Location:</b> select if the post-dispense should be purged in the source vessel of aspiration or should remain in the tips.</p> <p><b>Reuse Post-Dispense:</b> click on “Yes” if the post-dispense should remain in the tip and be used for the next aspiration.</p> <p><b>Dispense Type:</b> click on “Single” to aspirate and dispense the selected volume or on “Multi” to aspirate a multiple of the selected volume and perform several dispenses.</p> <p><b>Volume:</b> define whether the volume(s) to be dispensed (for every dispense step) is fix or variable. If variable, select each individual column and assign a volume. Set the pipetting direction. Additionally, to the dispense volume set also the <b>Pre-Dispense</b> and the <b>Post-Dispense</b> volume.</p> <p>Note: we recommend a pre- and post-dispense volume of 3-5 % of the nominal volume of the pipette.</p>
	Serial Dilution	<p><b>Last Aspiration:</b></p> <ul style="list-style-type: none"> <li>- “No”: the last aspiration is not performed; the liquid remains in the last row.</li> <li>- “Tip”: the last aspiration remains in the tips that are ejected into the Waste bin.</li> <li>- “Waste”: the last aspiration is purged in a dedicated liquid waste vessel. The location must be defined in “Pipetting location”.</li> </ul> <p><b>Volume:</b> the sample volume that will be transferred from well to well.</p>
	Multi Aspirate	<p><b>Aspirate Type:</b> click on “Single” to aspirate only one volume or on “Multi” to perform several aspirations consecutively.</p> <p><b>Volume:</b> define whether the volume to be aspirated (for every aspirate step) is fix or variable. If variable, select each individual column and assign a volume. Set the aspiration volume.</p>

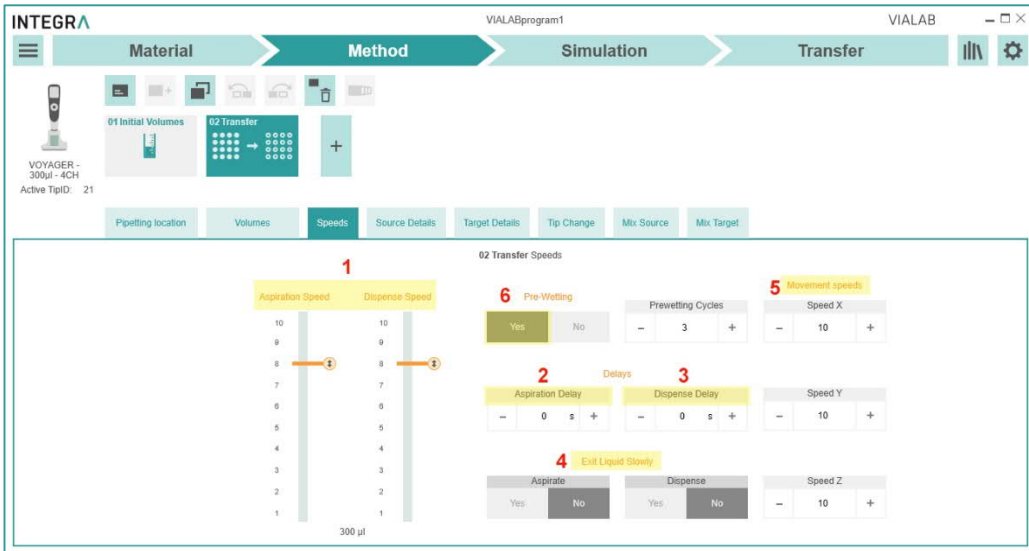


Parameter	Description
Speeds	<ul style="list-style-type: none"> <li>• <b>Aspiration and Dispense Speed (1)</b>: set the aspiration and dispense speeds for this step. Set the slider to the appropriate speed from 1 (slow) to 10 (fast).</li> <li>• If necessary, adapt parameters for viscous and volatile liquids:  <b>Aspiration Delay (2)</b>: a delay in seconds after aspiration before moving out of the liquid.  <b>Dispense Delay (3)</b>: a delay in seconds after dispensing before moving out of the liquid.  <b>Exit Liquid Slowly (4)</b>: the tips move slowly out of the liquid (for aspiration and for dispense).</li> </ul>
	 <p>The screenshot shows the INTEGRA VIALAB software interface. The top navigation bar includes 'Material', 'Method', 'Simulation', and 'Transfer'. The 'Method' tab is selected, showing a workflow with '01 Initial Volumes' and '02 Transfer'. The 'Speeds' sub-tab is active, displaying two vertical sliders for 'Aspiration Speed' and 'Dispense Speed', both set to 8. To the right, there are settings for 'Pre-Wetting' (6), 'Prewetting Cycles' (3), 'Movement speeds' (5) for Speed X (10), Speed Y (10), and Speed Z (10). Below these are 'Delays' (2 and 3) for Aspiration Delay and Dispense Delay, both set to 0 seconds. At the bottom, there are 'Exit Liquid Slowly' (4) settings for 'Aspirate' and 'Dispense', both set to 'No'.</p> <ul style="list-style-type: none"> <li>• <b>Movement Speeds (5)</b>: changes the speeds of the Tower in X-direction, of the Pipette arm in Z-direction or the Pipette holder in Y-direction [1 (slow) to 10 (fast)].</li> <li>• <b>Pre-Wetting (6)</b>: if pre-wetting of the tips is desired activate it and set the number of cycles. All fresh tips are pre-wetting with the following aspiration volume and speeds.</li> </ul>

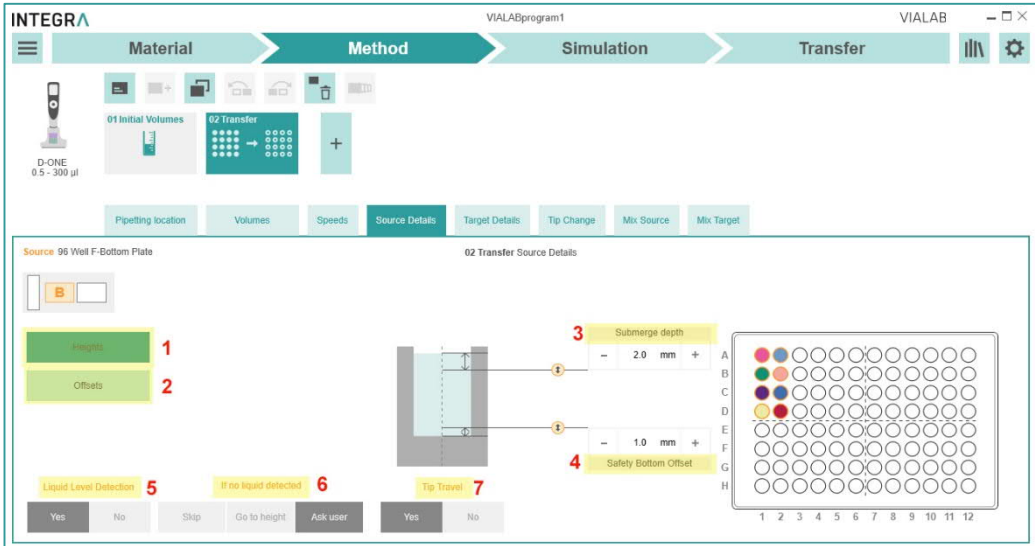
Parameter	Description
Source Details	The calculated liquid level in the well is displayed to ease height definition graphically.
Target Details	<p>Click on <b>Heights</b> (1) and move the line up or down to enter the heights for aspiration/dispense referring to the deck (deck height is 0 mm). The <b>Safety Bottom Offset</b> (4) defines the minimum distance of the tips to the bottom of the labware and is calculated automatically according to the labware dimensions and pipetting volume.</p> <p>If the <b>Heights</b> are variable, define the pipetting heights for every column.</p> <p><b>Tip Travel</b> (7) is the movement of tip during one aspiration/dispense to compensate the decrease/increase of the liquid level. It is recommended for high pipetting volumes in vessels with small diameter. If activated, aspiration/dispense starts at the <b>Start Height / Submerge Depth</b> (3) and ends at the <b>End Height / Safety Bottom Offset</b> (4), which is calculated automatically according to the labware dimensions and pipetting volume. If Tip Travel is deactivated the aspiration/dispense is performed at a constant <b>Height</b> (1).</p>
	 <p><b>D-ONE only:</b></p> <p><b>Liquid Level Detection</b> (5) is on per default and the <b>Submerge depth</b> (3) is displayed on top of the well. The liquid level is detected by the D-ONE automatically. If "No" is selected, the pipetting start height must be set manually.</p> <p>Define the error handling, <b>If no liquid detected</b> (6). "Skip" automatically continues with the next source well. "Go to height" continues by aspirating from a pre-defined pipetting height. "Ask user" stops and asks the user to select one of these options: Skip, Repeat and Go to height.</p> <p><b>Note:</b> it is possible to fine-tune the Pipetting Height / Start Height directly on the pipette to optimize the heights according to the used liquid (see ASSIST PLUS Operating Instructions chapter 5.4 "Adjust heights and positions").</p>





**Note:** LLD feature exclusive for the D-ONE modules and not for VIAFLO/VOYAGER pipettes.

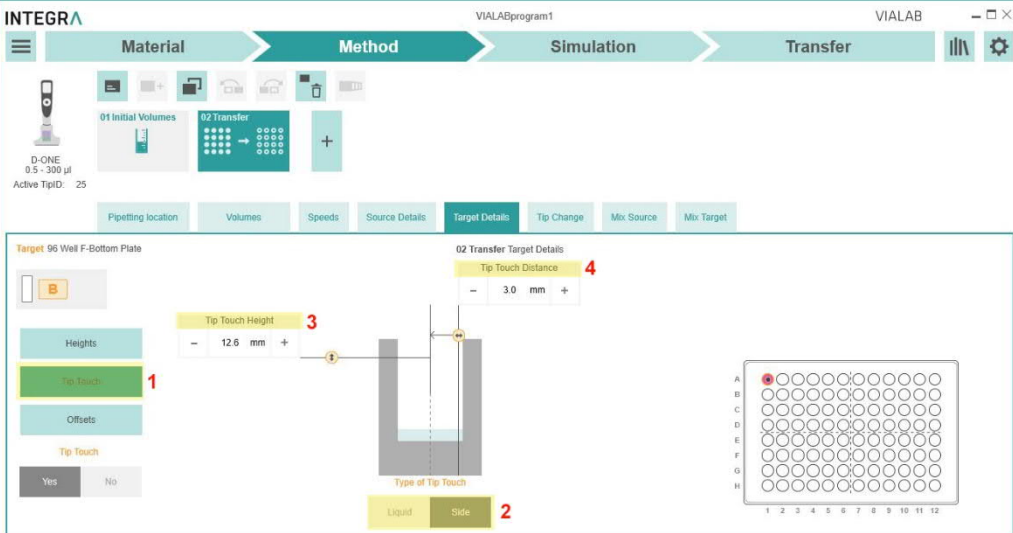Parameter	Description
Tip Touch	<p>Set a tip touch (1) after dispensing to remove drops that may cling on the pipette tips. Set whether the tip touch should be performed into the <b>Liquid</b> or on the <b>Side</b> (2) of the wells. Enter the <b>Tip Touch Height</b> (3) and the <b>Distance</b> (4) to the vessel wall.</p> 

**Note:** If red exclamation marks appear, the program contains errors. Click on “Validation messages” for further details.

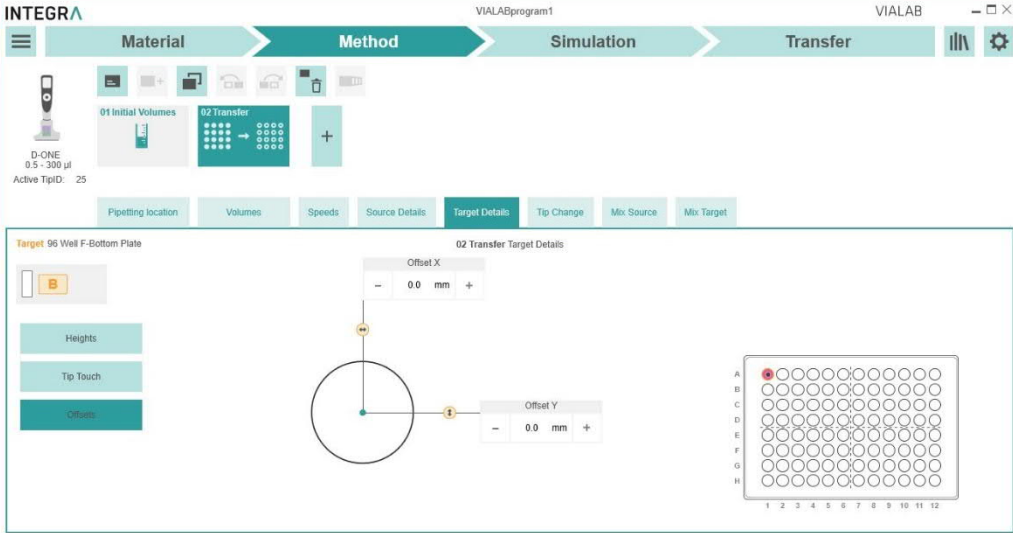Parameter	Description
Offset	 <p>Define a well <b>Offset</b>, if you want to pipette off the well center or even out of wells, this is helpful for processing e. g. for transwell plates.</p>

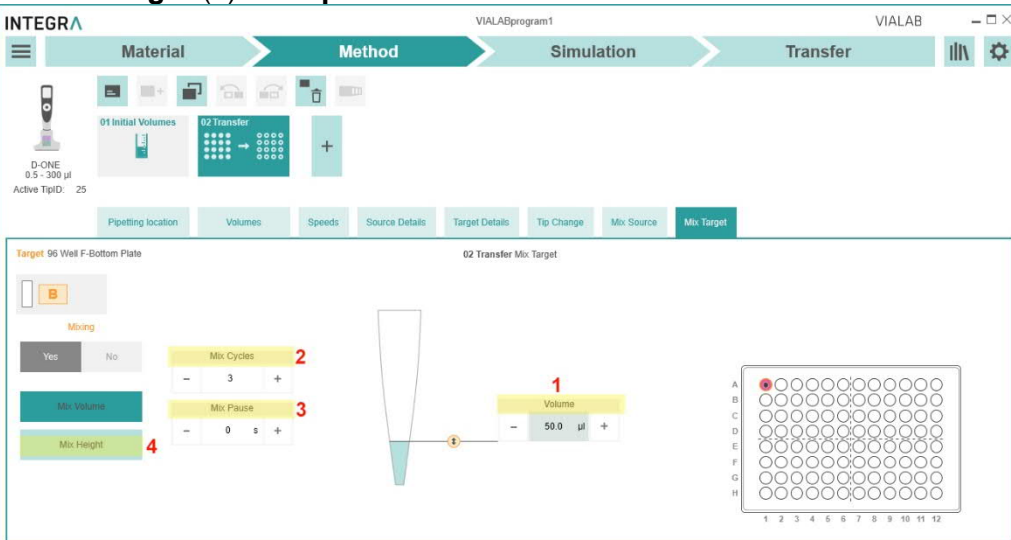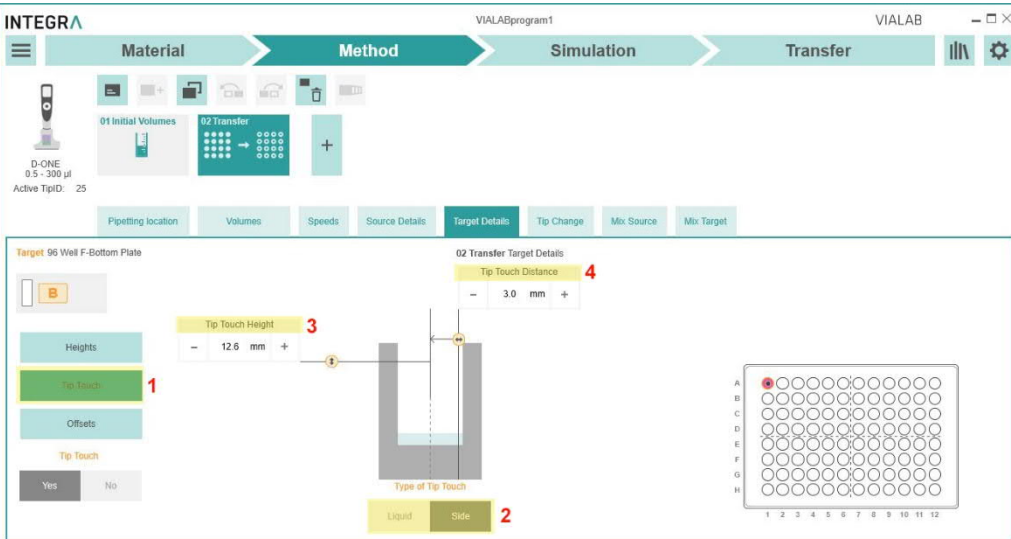


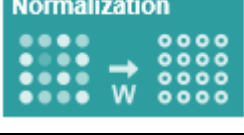

Parameter	Description
Tip Change	<p>Set when tips should be changed. The Tip change options can vary from step to step in VIALAB. Here some examples:</p> <p><b>Continue next step without change:</b> same tips will be used for next step.  <b>After step complete:</b> change tips after completing the current step  <b>After each transfer/dispense/dilution:</b> change tips after each dispense  <b>Before each new aspiration:</b> change tips before each new aspiration</p> <p>If the “Tip Ejection Type” is activated in the Program Settings (see chapter 0), the tips can be either positioned on the ejector ledge before ejection (<b>Ledge</b>) or are directly ejected into the middle of the waste bin (<b>Waste</b>). Tips may stack up and the waste bin needs to be emptied more often.</p>
Mix	<p>Define, if the sample should be mixed in the source well before aspiration and in the target well after dispensing. Define the mix <b>Volume</b> (1), the number of <b>Cycles</b> (2), if a <b>Pause</b> (3) should be performed in between mixing cycles. Define the <b>Mix Height</b> (4) and <b>Speed</b>.</p> 
Tip Touch	<p>Set a tip touch (1) after dispensing to remove drops that may cling on the pipette tips. Set whether the tip touch should be performed into the <b>Liquid</b> or on the <b>Side</b> (2) of the wells. Enter the <b>Tip Touch Height</b> (3) and the <b>Distance</b> (4) to the vessel wall.</p> 



**Note:** If red exclamation marks appear, the program contains errors. Click on “Validation messages” for further details.

## 6.4 Method steps only for D-ONE pipetting module

These additional steps are only available for D-ONE pipetting module (PN 4531, 4532):

Symbol	Step	Description
	Master Mix	Add multiple master mix reagent components from individual tubes to one target tube or well. A Master Mix reaction calculator generates reagent volumes required.
	Hit Picking	Individual wells from multiple source plates are transferred to one target plate (condense the hits). Create and import a worklist containing the source plate and sample location.
	Normalization	Sample of different concentrations are diluted to a final concentration. Pipetting location, sample identifier and concentration data are entered manually by the user.
	Normalization Worklist	Sample of different concentrations are diluted to a final concentration. Pipetting location, sample identifier and concentration data is imported from a worklist file.
	Worklist	Pipetting steps are created based on an imported worklist. Allows higher pipetting flexibility useful for e.g., NGS Library Pooling and any sample transfers imported in a table format.

### 6.4.1 Importing worklists

D-ONE module allows importing worklists that either follow a specific worklist format or have been created using the worklist templates provided by INTEGRA. The pipetting steps Hit Picking, Normalization Worklist and Worklist allow import of worklists with max. 384 lines.

#### 6.4.1.1 Worklist format

Hit Picking worklist:

File format:	.csv
Column separator:	semicolon
Number of columns	6
Column header	SampleID, SourcePlateID, SourceWell, TargetPlateID, TargetWell, TransferVolume

Find below an example of a Hit Picking worklist:

```
SampleID;SourcePlateID;SourceWell;TargetPlateID;TargetWell;TransferVolume
Sample 01;Source 1;C5;TARGET;A1;25
Sample 02;Source 1;H11;TARGET;B1;25
Sample 03;Source 1;D3;TARGET;C1;25
Sample 04;Source 1;C8;TARGET;D1;25
Sample 05;Source 1;A7;TARGET;E1;25
```

Normalization worklist:

File format:	.csv
Column separator:	semicolon
Number of columns	6
Column header	SampleID, SourcePlateID, SourceWell, TargetPlateID, TargetWell, InitialConcentration

Find below an example of a Normalization worklist:

```
SampleID;SourcePlateID;SourceWell;TargetPlateID;TargetWell;InitialConcentration;
01 DNA sample;Sample Plate;A1;Destination Plate;A1;39;
02 DNA sample;Sample Plate;B1;Destination Plate;B1;38;
03 DNA sample;Sample Plate;C1;Destination Plate;C1;37;
04 DNA sample;Sample Plate;D1;Destination Plate;D1;77;
05 DNA sample;Sample Plate;E1;Destination Plate;E1;89;
```

Worklist:

File format:	.csv
Column separator:	semicolon
Number of columns	6
Column header	SampleID, SourceDeckPosition, SourceWell, TargetDeckPosition, TargetWell, TransferVolume [µl]


Find below an example of a Worklist:

```
SampleID;SourceDeckPosition;SourceWell;TargetDeckPosition;TargetWell;TransferVolume [µl]
Sample 1;B;A1;C;A1;50
Sample 2;B;B1;C;B1;100
Sample 3;B;C1;C;C1;125
Sample 4;B;D1;C;D1;25
Sample 5;B;E1;C;E1;100
Sample 6;B;F1;C;F1;50
Sample 7;B;G1;C;G1;75
```

#### 6.4.1.2 Create a worklist with INTEGRA's template


- If the worklist is already defined, just press **Import (1)** to import the worklist into VIALAB.
- To create a new worklist, press on **Templates (2)** and you are automatically directed to the VIALAB folder "MyWorklists" containing templates and sample work lists for Hitpicking and Normalization.



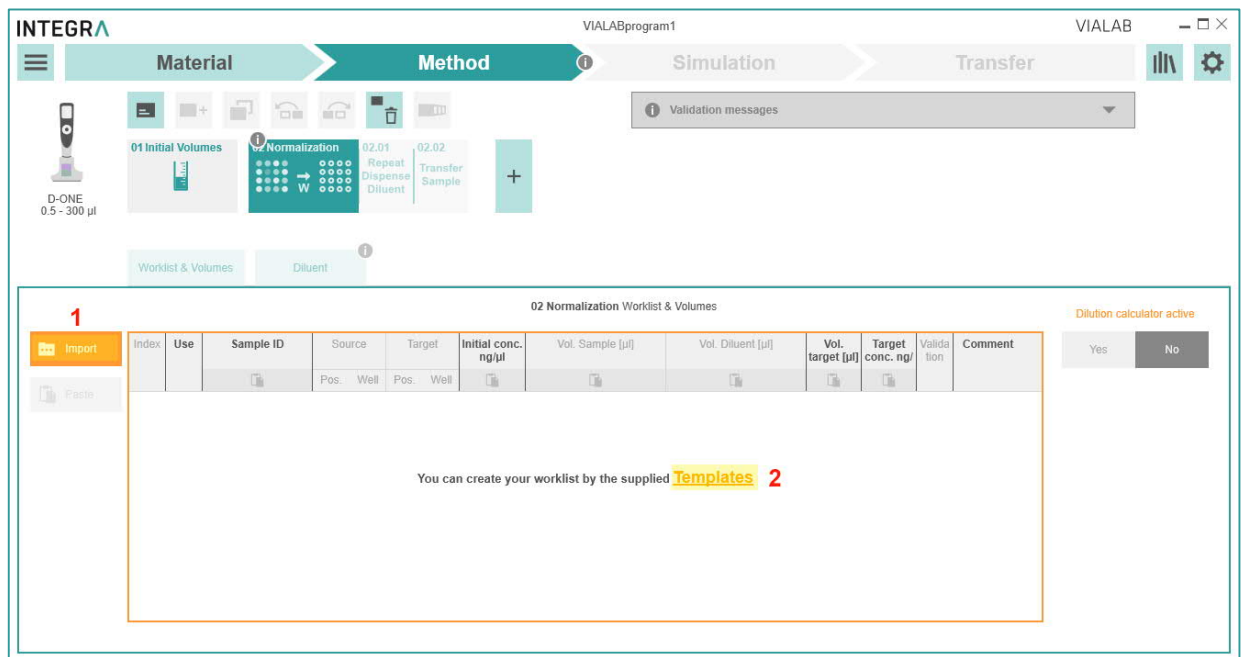


**Note:** Always use the templates provided by INTEGRA to generate a worklist.

- Open the corresponding template, e.g. the normalization template:

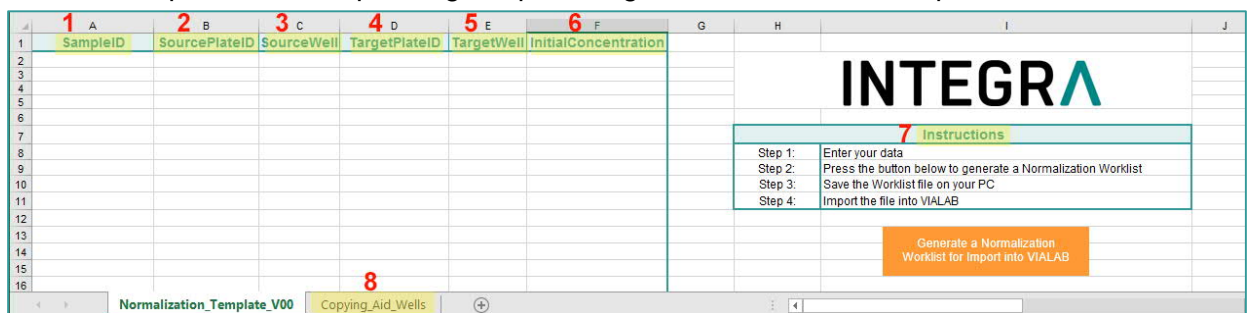

- 1 SampleID:** name of sample, e.g. barcode.
- 2 SourcePlateID:** name of source plate, e.g. barcode or source plate 1, source plate 2, etc.
- 3 SourceWell:** position of source well, e.g. A1, B3, etc.
- 4 TargetPlateID:** name of target plate, e.g. target plate.
- 5 TargetWell:** position of target well.
- 6 TransferVolume [ul]:** transfer volume.
- 7 Instructions:** follow the instructions to generate the Normalization template.
- 8 Copying\_Aid\_Wells:** list of well numbers which can be used to copy source and target wells into the template.

- Copy data from e.g. another excel file column by column into the template file, or use Copying\_Aid\_Wells (8):

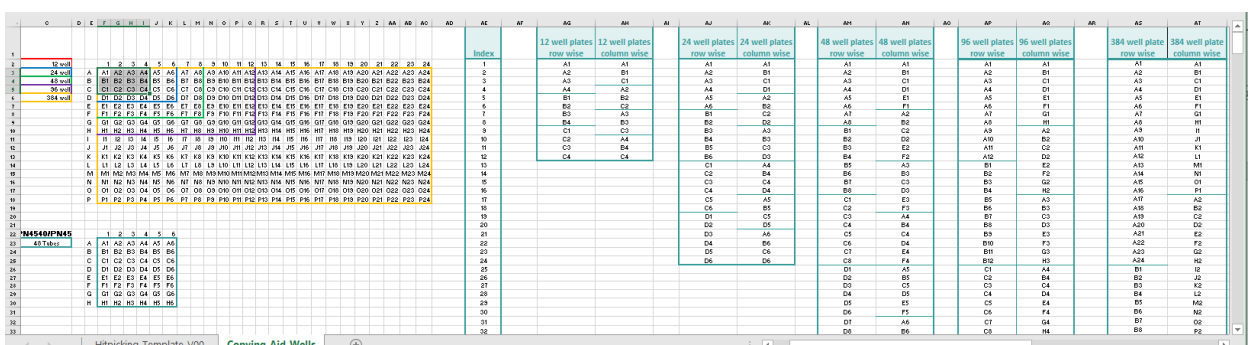


- After the template is filled out press on the orange button "Generate a Normalization Worklist", and a worklist.csv file is automatically generated. Name the file and save on e.g. your local drive.
- Go back to VIALAB and import the previous defined worklist.csv file.




## 6.4.2 Hit Picking

- Start your method by adding a **Hit Picking (1)** step.

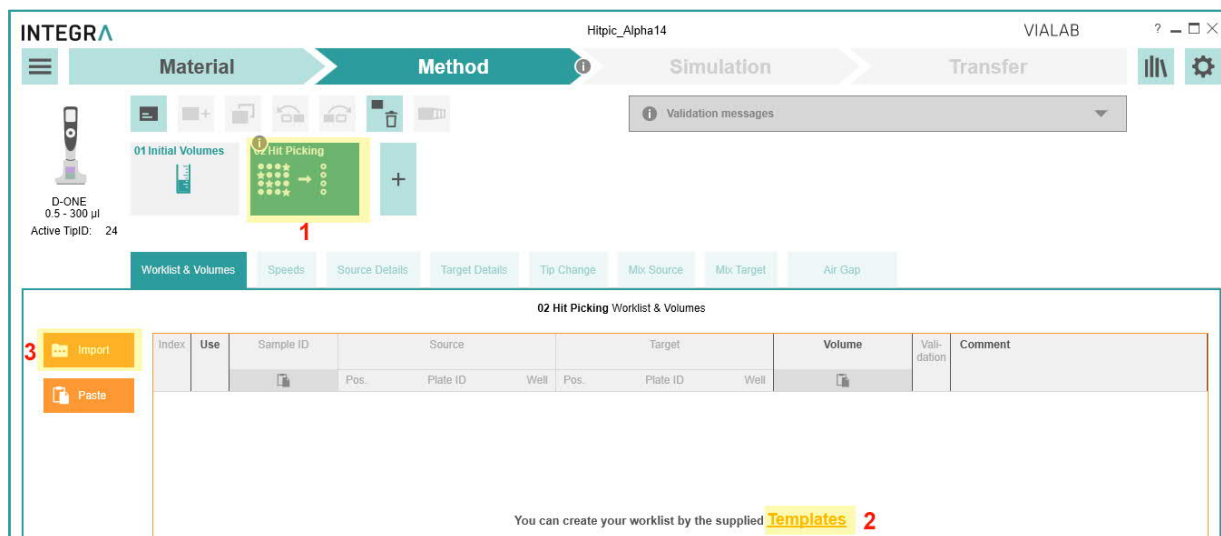


- Click on **Templates (2)** to open an EXCEL file and define a worklist with source and target wells, see 6.4.1. After pressing the orange button in the template, a worklist in a csv file format necessary for import into VIALAB is generated. Click on **Import (3)** to load the worklist.
- You are requested to assign the deck position where the target plate is placed.

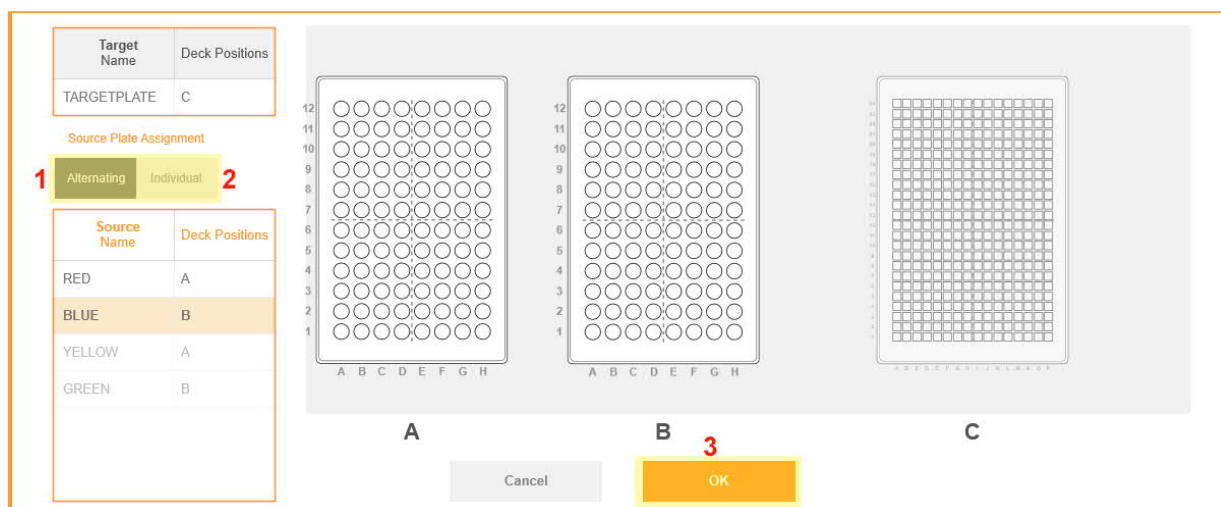

If the source plates should be processed in a fix pattern, select **Alternating (1)** and just assign plates and corresponding deck positions of the first layout. Click on **Individual (2)** if the source plates should be assigned individually to each deck position. Press **OK (3)**.


- The worklist is displayed:

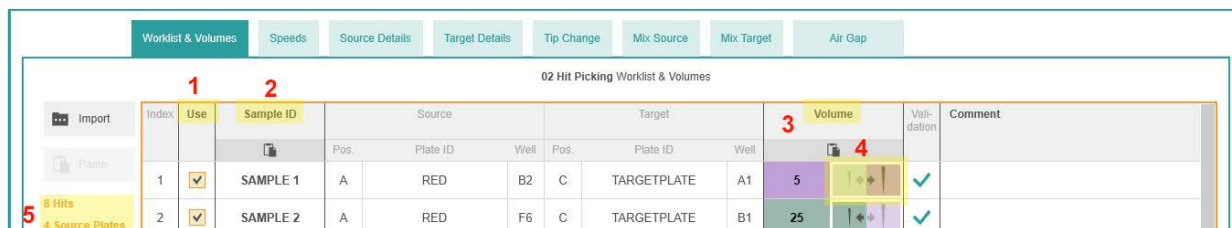

By removing the checkmark under **Use (1)**, you can exclude samples from processing. You can also change the sample ID (2), the volume (3) and the tip type (4) and. The number of hits and source plates are shown at the side (5).

- Adjust the other parameters like Speeds, Source and Target Details, Tip Change.
- Under **Source Details** (1) **Liquid Level Detection** (LLD, 2) and **Tip Travel** (4) are default switched on. Consequently, the aspiration height and the initial volume must not be defined.

**Note:** LLD is only possible with new and dry tips.

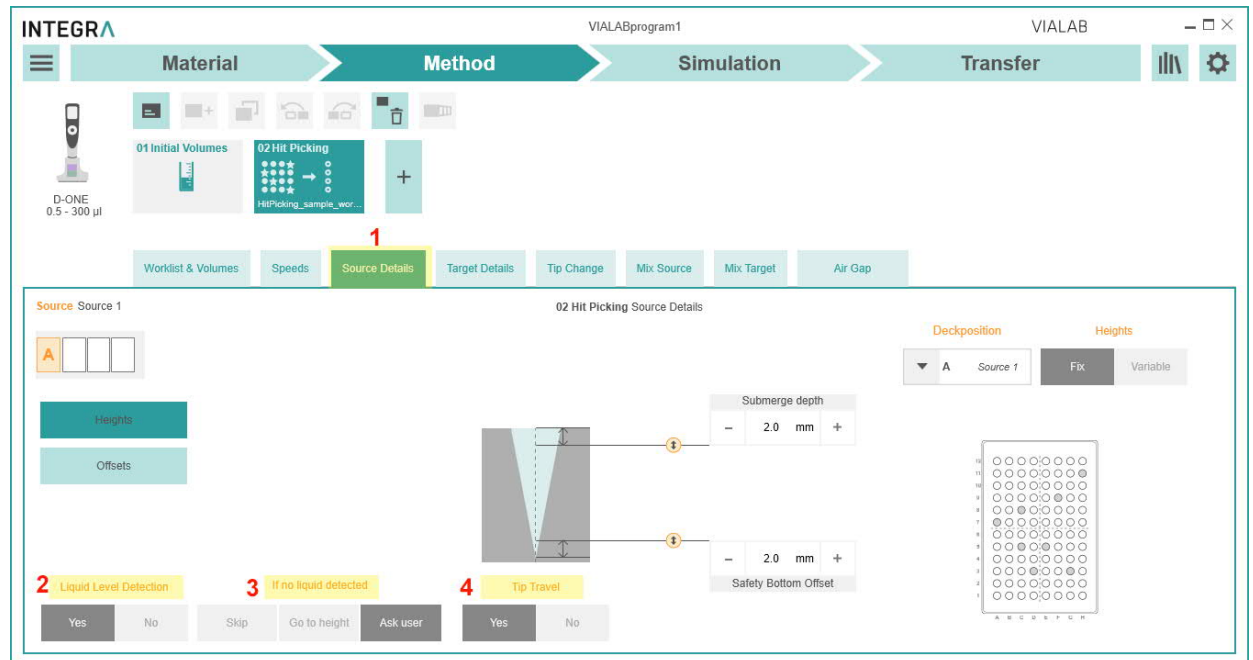


Define the error handling, **If no liquid detected** (3), as described under Method, see 6.3.

### 6.4.3 Normalization

With this step the samples which need to be normalized must be chosen in the graphics.

- Add a **Normalization** (1) step to your method. Define the source (2), target (3) wells and then the diluent (4).

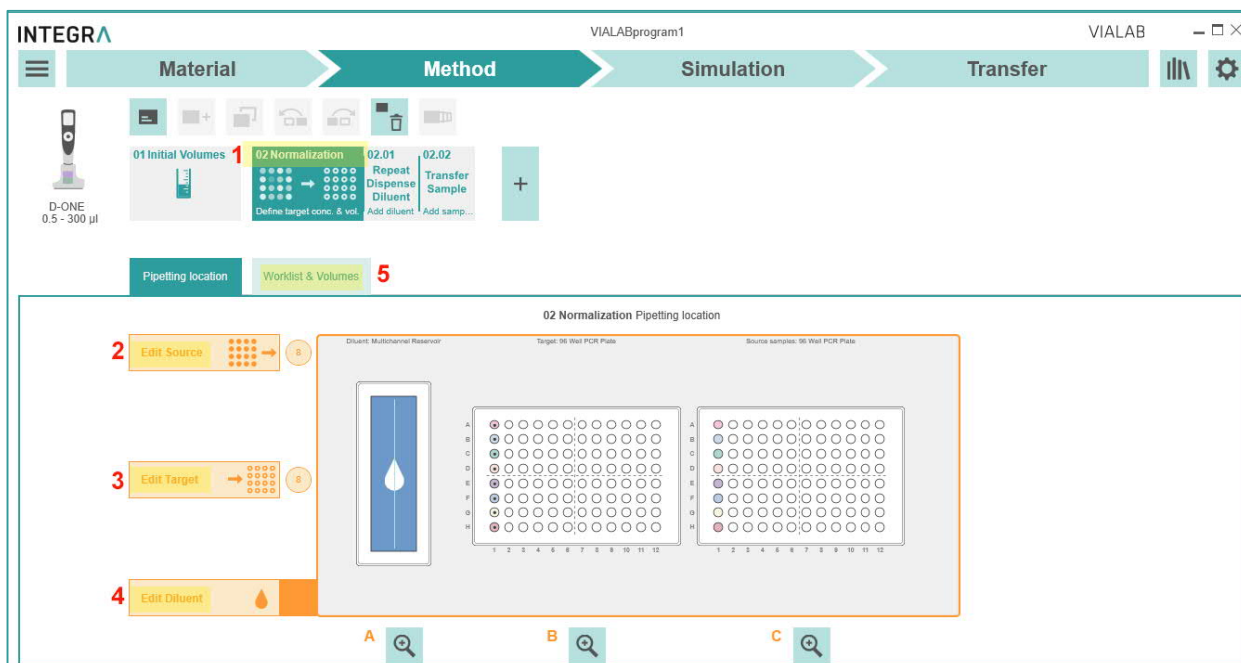

- Under **Worklist & Volumes** (5) the worklist is displayed:

The screenshot shows the 'Worklist & Volumes' section of the INTEGRA VIALAB software. It displays a table with 6 rows of data for DNA samples. The table has columns for Index, Use, Sample ID, Source, Target, Initial conc. ng/μl, Vol. Sample [μl], Vol. Diluent [μl], Vol. target [μl], Target conc., Validation, and Comment. A 'Dilution calculator active' button is visible on the right side of the table.

Index	Use	Sample ID	Source	Target	Initial conc. ng/μl	Vol. Sample [μl]	Vol. Diluent [μl]	Vol. target [μl]	Target conc.	Validation	Comment
1	<input checked="" type="checkbox"/>	01 DNA sample	C A1	B A1	39	20.5	59.5	80	10	✓	
2	<input checked="" type="checkbox"/>	02 DNA sample	C B1	B B1	38	21	59	80	10	✓	
3	<input checked="" type="checkbox"/>	03 DNA sample	C C1	B C1	37	21.5	58.5	80	10	✓	
4	<input checked="" type="checkbox"/>	04 DNA sample	C D1	B D1	77	10.39	69.5	79.89	10	✓	
5	<input checked="" type="checkbox"/>	05 DNA sample	C E1	B E1	89	8.99	71	79.99	10	✓	
6	<input checked="" type="checkbox"/>	06 DNA sample	C F1	B F1	48	16.67	63.5	80.17	10	✓	

- By removing the checkmark under Use (1), you can exclude samples from processing. Copy or type sample details, i.e. the name/barcode (Sample ID, 2) and the initial concentration (3).
- If you switch on the **Dilution calculator** (4), you can adjust the target volume and concentration. Click on **Calculate all samples** (5) to automatically calculate all parameters (volumes, pipette tips) automatically.

- The step **Repeat Dispense Diluent (1)** first transfers the diluent, followed by **Transfer Sample (2)**. Both steps allow to set pipetting parameters such as Speeds.

**INTEGRA** VIALABprogram1 VIALAB

**Material** **Method** **Simulation** **Transfer**

D-ONE 0.5 - 300 µl

01 Initial Volumes 02 Normalization 02.01 Repeat Dispense Diluent 02.02 Transfer Sample

Define target conc. & vol. Add diluent Add samp...

**1** **2**

**Speeds** Volumes Source Details Target Details Tip Change Mix Source

**02.01 Repeat Dispense Diluent Speeds**

**Aspiration Speed** **Dispense Speed**

10 9 8 7 6 5 4 3 2 1

10 9 8 7 6 5 4 3 2 1

300 µl

**Pre-Wetting**

Yes No

**Delays**

**Aspiration Delay** **Dispense Delay**

- 0 s + - 0 s +

**Exit Liquid Slowly**

**Aspirate** **Dispense**

Yes No Yes No

**Movement speeds**

**Speed X**

- 10 +

**Speed Y**

- 10 +


**Speed Z**

- 10 +


#### 6.4.4 Normalization Worklist


With this step the samples which need to be normalized must be previously defined in a worklist and imported.

- Add a **Normalization Worklist** (1) step to your method. Click on **Import** (2) to load the worklist file, see 6.4.1.

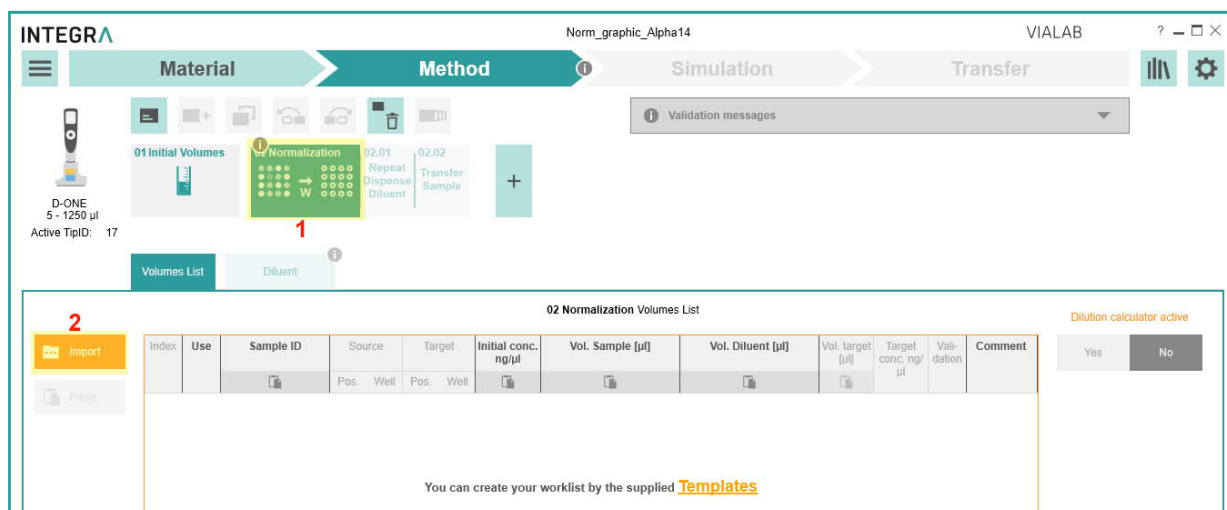

- Assign where the target and source plates are placed on deck. Then assign the position of the diluent.

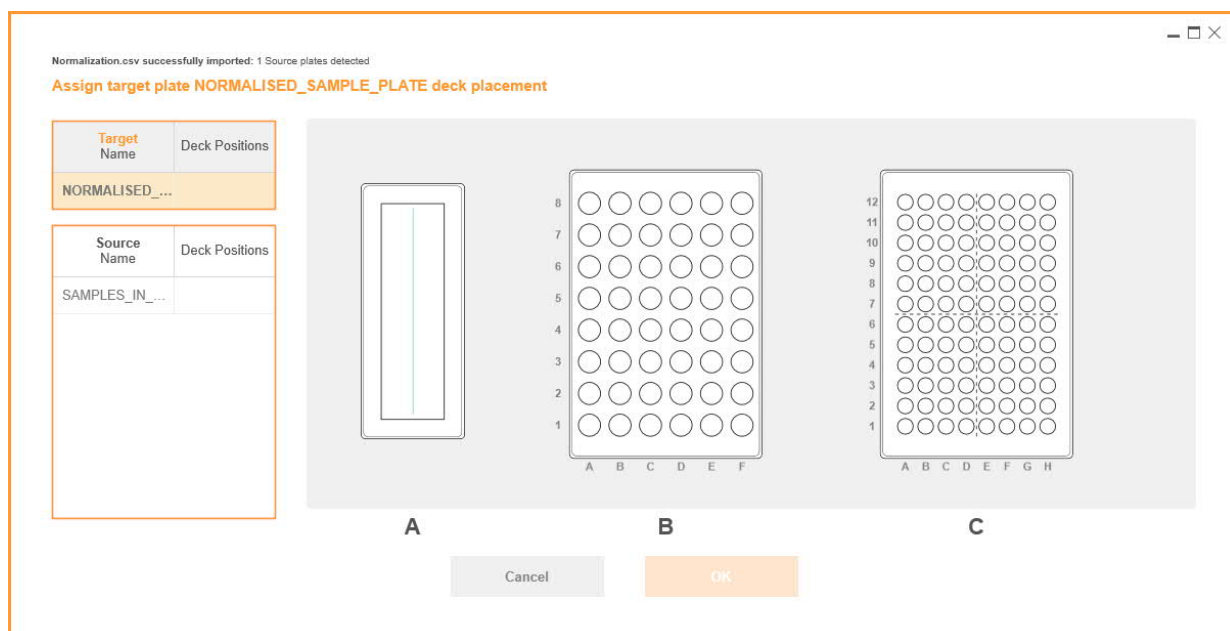

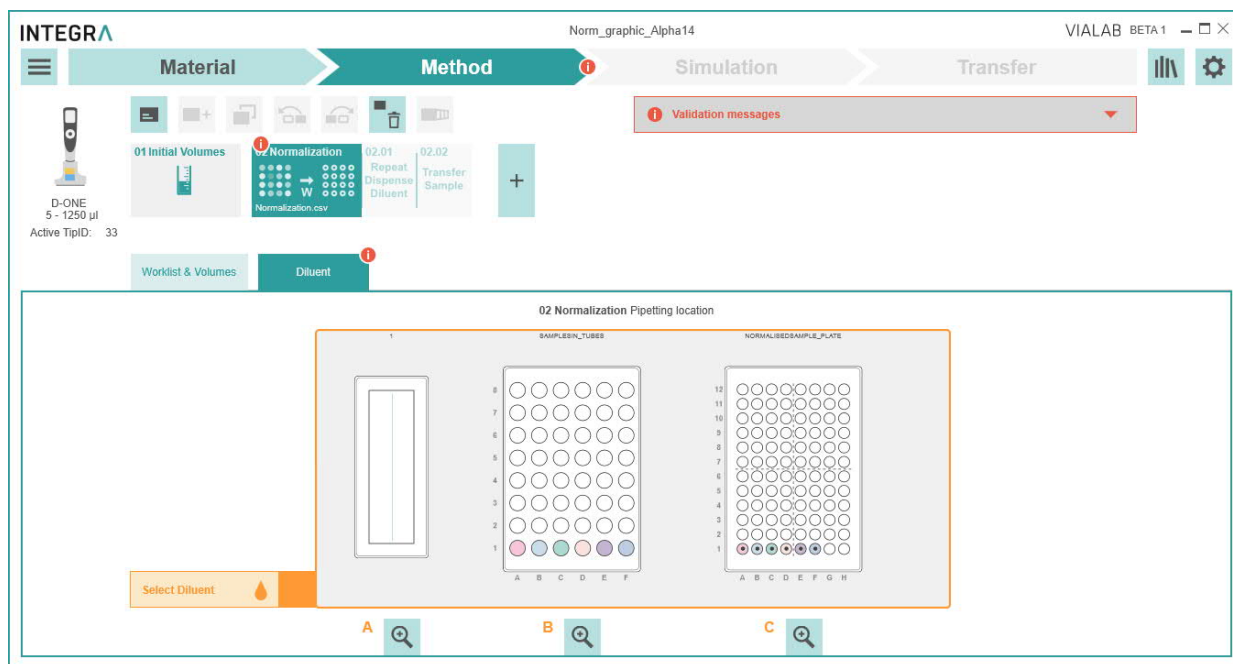

Click on the diluent, i.e. the reservoir in the example above.

- Under **Volumes List** (1) the worklist is displayed. By removing the checkmark under Use (2), you can exclude samples from processing.

The screenshot shows the 'Volumes List' table and the 'Dilution calculator' panel. The table lists 6 samples with their respective volumes and concentrations. The 'Dilution calculator' panel is active, showing target volume and concentration settings.

Index	Use	Sample ID	Source	Target	Initial conc. ng/μl	Vol. Sample [μl]	Vol. Diluent [μl]	Vol. target [μl]	Target conc. ng/μl	Validation	Comment
1	<input checked="" type="checkbox"/>	SAMPLE 1	B A1	C A1	86	20	180	200	8.6	✓	
2	<input checked="" type="checkbox"/>	SAMPLE 2	B B1	C B1	50	34.4	166	200.4	8.6	✓	
3	<input checked="" type="checkbox"/>	SAMPLE 3	B C1	C C1	30	57.3	143	200.3	8.6	✓	
4	<input checked="" type="checkbox"/>	SAMPLE 4	B D1	C D1	55	31.3	169	200.3	8.6	✓	
5	<input checked="" type="checkbox"/>	SAMPLE 5	B E1	C E1	63	27.3	173	200.3	8.6	✓	
6	<input checked="" type="checkbox"/>	SAMPLE 6	B F1	C F1	52	33.1	167	200.1	8.6	✓	

**3 Dilution calculator active**

Yes No

Target volume

Fixed Variable

20 Sample Volume (μl)

200 Target volume (μl)

Target concentration


8.6 ng/μl

**4 Calculate all samples**


If you switch on the **Dilution calculator** (3), you can adjust the target volume and concentration. Click on **Calculate all samples** (4) to automatically calculate all parameters (volumes, pipette tips) automatically.

## 6.4.5 Master Mix

- Define the location of the source (1) and the target (2) tubes.

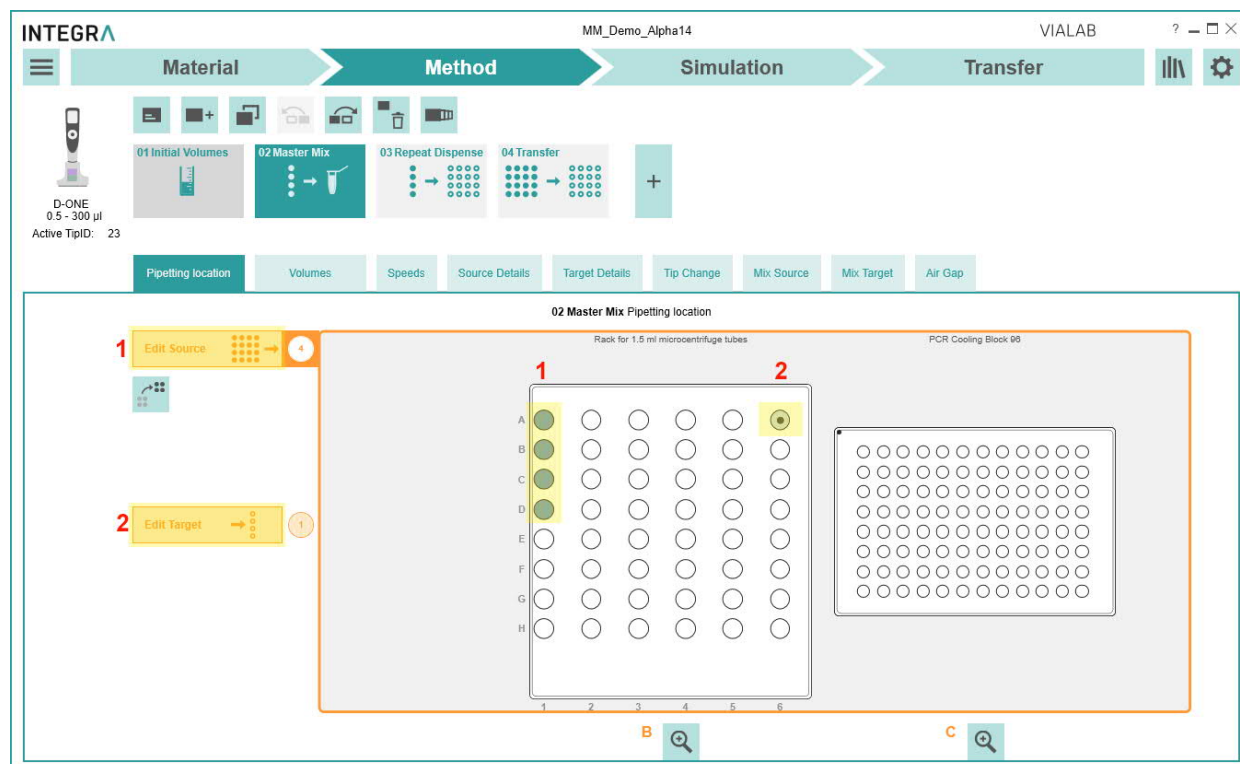

- Under **Volumes** name the reagents (1), e.g. water, dNTP, enzyme, primer, and enter the Volume (2) for a single reaction.

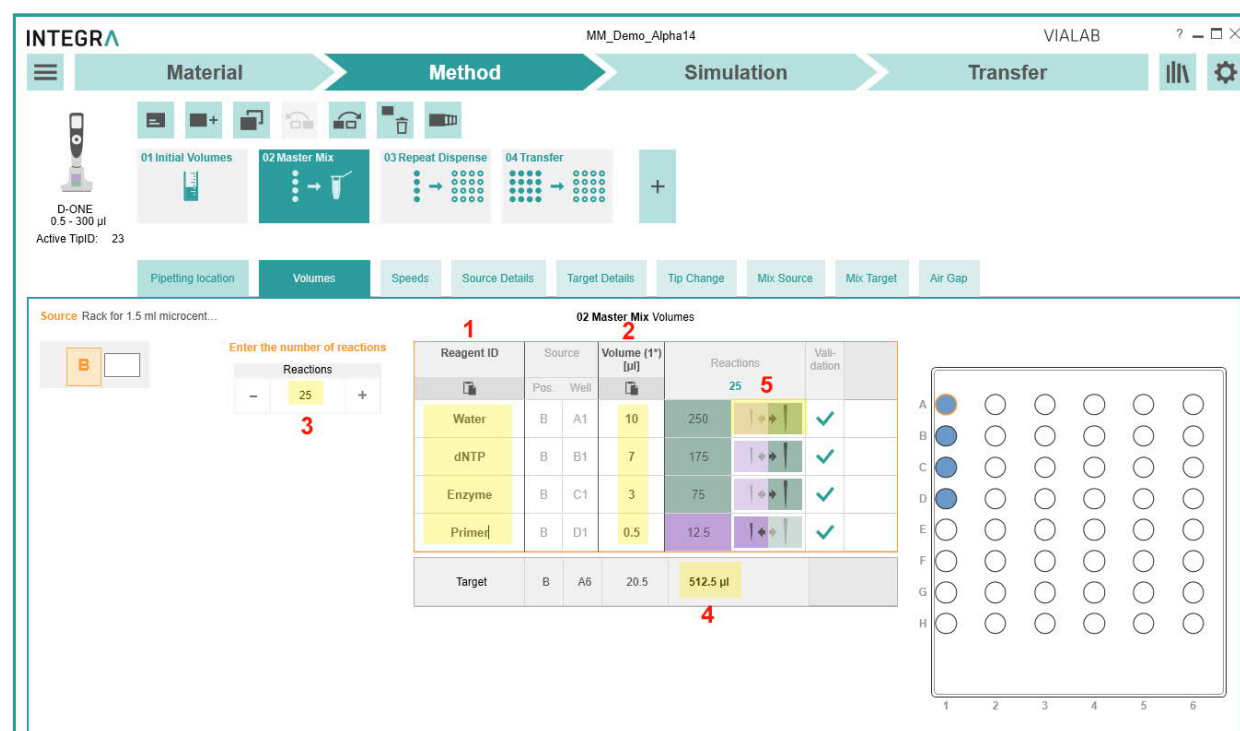

Enter the number of reactions (3). The total reagent volumes and the Target (4) volume are calculated automatically. The correct tip size is also selected automatically by VIALAB. If desired, the tip type can be manually changed by the user (5).



**Note:** If the calculated volume exceeds the tip volume, the volume is dispensed multiple times with the same tip, indicated by a multiplier (\*, **1**).

Pipetting location Volumes Speeds Source Details Target Details Tip Change Mix Source Mix Target

Source Rack for 1.5 ml microcent...

02 Master Mix Volumes

Enter the number of reactions

Reactions: - 48 +

Reagent ID	Source	Volume (1*) [µl]	Reactions	Validation
Water	B A1	10	480 (240 µl*2)	✓
dNTP	B B1	7	336 (168 µl*2)	✓
Enzyme	B C1	3	144	✓
Primer	B D1	0.5	24	✓
Target	B A6	20.5	984 µl	

## 6.4.6 Worklist

Use the Worklist step to set up a customized program based on an imported worklist, see 6.4.1.

INTEGRA Hitpic\_Alpha14 VIALAB

Material Method Simulation Transfer

D-ONE 0.5 - 300 µl Active TipID: 24

01 Initial Volumes Worklist

Validation messages

Worklist & Volumes Speeds Source Details Target Details Tip Change Mix Source Mix Target Air Gap

02 Worklist Worklist & Volumes

Import Paste Overview

Index	Use	Sample ID	Source	Target	Volume	Validation	Comment
			Pos. Well	Pos. Well			

You can create your worklist by the supplied [Templates](#)




## 6.5 Simulation

In the simulation window the deck with the defined labware is displayed. On the top all steps are listed. Click on the play button to start the simulation of your program. Watch the whole simulation or navigate between steps or stop the simulation (1).

To speed up the simulation to 2-, 4- or 8-fold speed, select the appropriate button (2). An estimation of the Total Time and the Total Tip Consumption is listed (3).

The simulation is displayed as top (4) and side (5) view of the deck.

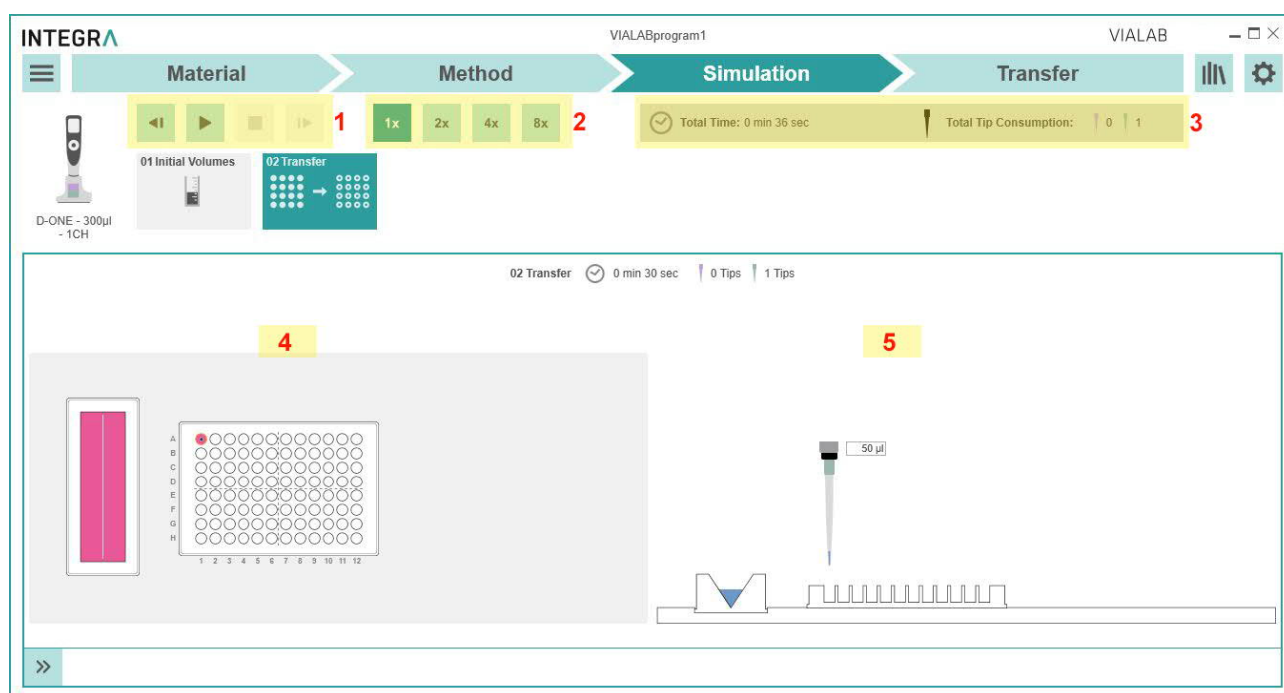

For all error-free programs a report is generated by clicking on the Main Menu and selecting "Report". A report in pdf format with all settings of the program is generated.

## 6.6 Transfer

Transfer your VIALAB program to your pipette and run it from your pipette. Alternatively connect the ASSIST PLUS via USB to the PC, connect the Pipette and start a program from VIALAB directly.

### 6.6.1 Connect pipette or pipetting module to the PC

#### Connection via charging/communication stand (#4211)

Connect the USB cable to the stand and a free USB port on your computer. It is recommended to connect the stand to the mains power supply. Then hang the pipette or the D-ONE module onto the stand.

Go to the “Toolbox” of the pipette/D-ONE. In the menu “Communications” select “USB”. Stay in this menu for transfer.

#### Connection via communication module for pipette (#4221)


Pipette or D-ONE and PC need to be paired for the first connection. Please refer to MS Windows Help to learn how to pair Bluetooth devices. The following pairing code may be required: 12345.

Go to the “Toolbox” of the pipette. In the menu “Communications” select “via ASSIST PLUS”.

### 6.6.2 Transfer

In VIALAB switch to the section “Transfer” for saving the program to the pipette.

Upon establishing the connection successfully, the pipette details, such as pipette model, serial number, and firmware version, are displayed on the right side. All VIALAB programs already saved on the pipette are also listed.

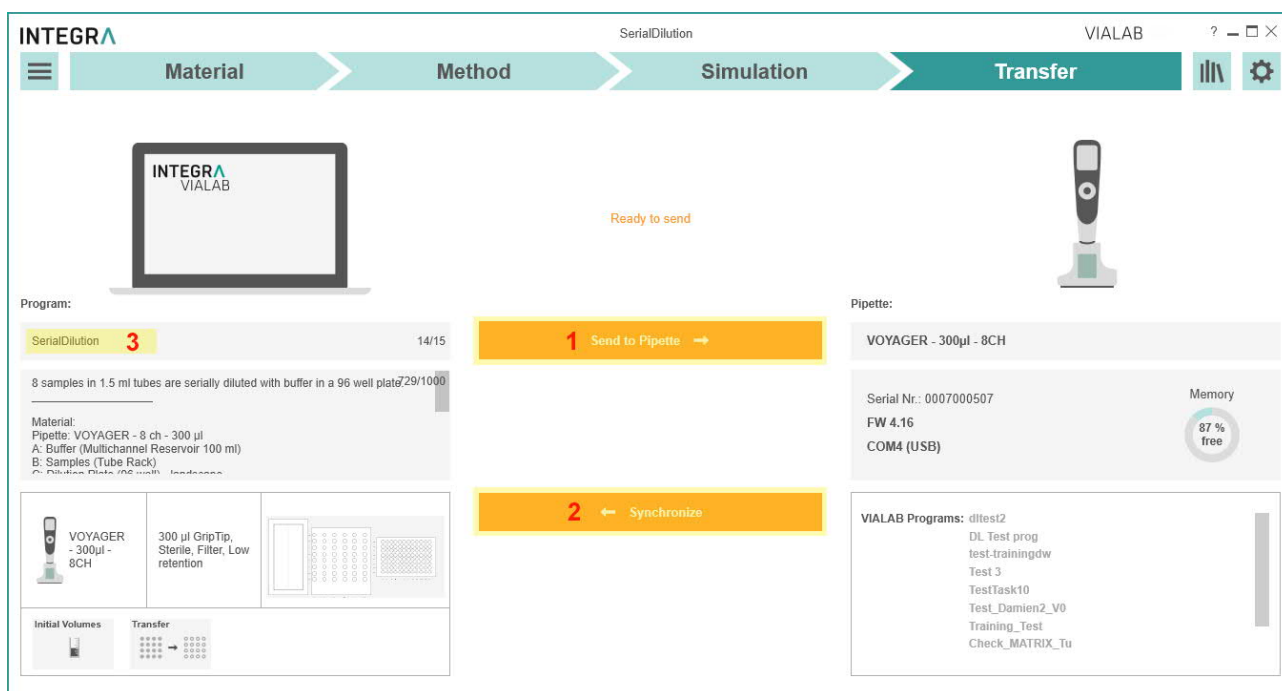

Send to Pipette (1) starts the transfer of the current program to the pipette. If certain parameters of the program are directly modified on the pipette (offsets, pipetting heights), synchronization of the program to VIALAB is possible (2).


To change the program display name on the pipette type in an alternative program name with max. 15 characters (3).

Successful saving or synchronization of the program to the pipette / to VIALAB is indicated by a temporary message at the bottom of the window.

## 6.7 Run a program directly from VIALAB

To run a program directly from VIALAB, connect your ASSIST PLUS to your PC via USB. Then connect your pipette via ComModule to the ASSIST PLUS.

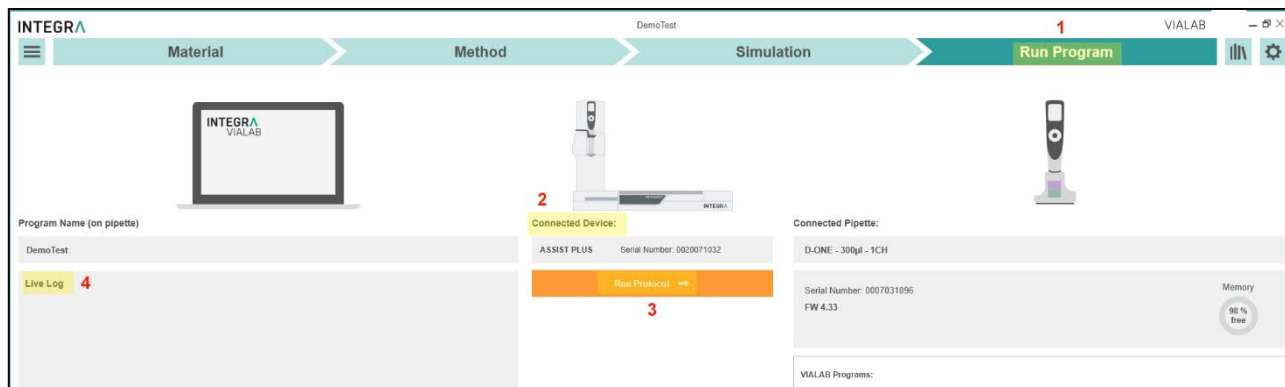


Enable a connection on the pipette under Toolbox – Communication, select PC Control A+.



As soon as the connection is successfully established the Transfer tab changes to Run Program (1) and the connected device is displayed with serial number (2).

Press Run Protocol (3). The VIALAB program is transferred automatically to the pipette. Under Live Log (4) the actual steps are displayed.

After the Run is done, the run report is generated automatically, showing all information such as serial numbers of ASSIST PLUS, pipette, user information, deck layout, method.

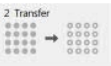VIALABprogram1 Report		INTEGRA VIALAB	
Program Name		VIALABprogram1.iaa	
Program Name (on pipette)		Test	
Last Saved Date:		24. Mar 2022	
Last Save Operator:			
Instrument - Serial Number		0020081124	
Pipette - Serial Number		0007031096	
Tip Type (PN 6434) Lot Nr.:			
Run Operator:		jan-peter.baldin	
Run Date:		24.Mar.2022	
Run Start Time:		11:30	
Run End Time:		11:31	
Notes:		11:30:57 : Run started 11:31:00 : Transfer (Step 02) 11:31:26 : Run finished	
Signature:			
Overview Method			
<div><p>D-ONE - 300µl - 1CH</p></div> <div><p>1 Initial Volumes</p></div> <div><p>2 Transfer</p></div>			
Total Time:		0 min 29 sec	
Total Tip Consumption:		Tip 300µl: 1	
Deck Layout			

## 7 Labware library

You can access the labware library by clicking on the Library icon in the top right corner.

Within the labware library, you can define the dimensions of decks, reservoirs, plates, tube racks and COMBI Bases/Systems. Pipettes, GripTips and reservoirs from INTEGRA are already defined and cannot be modified.

### 7.1 Reservoirs

In the reservoir menu all predefined plates are listed. Select a desired property, e.g. Type, to quickly show already defined reservoirs of this property.

Reservoir	Name	Volume	Manufacturer	Part Number	Description	Favorites
	12 Column Polypropylene Reservoir	21 ml	INTEGRA	6361, 6362		★
	8 Row Polypropylene Reservoir	32 ml	INTEGRA	6371, 6372		★

- 1 Reservoir:** the icon shows the number and shape of the wells and well bottom.
- 2 Name:** indicates the number and shape of the wells and well bottom, as well as other properties such as half-area.
- 3 Volume:** shows the maximal well volume of the reservoir or compartment as defined by the manufacturer.
- 4 Manufacturer:** name of the reservoir manufacturer
- 5 Part Number:** part number in accordance to the reservoir manufacturer
- 6 Description:** optional
- 7 Favorites:** select to mark as favorite. The material section in the program set-up has a filter to show only favorites.
- 8 Tools:** Edit/View details, Copy, Export or Delete the selected labware. It is not possible to edit or delete INTEGRA items. To modify other predefined labware, create a copy first and then edit the copy.

To add new reservoirs to the library, click on “Define new item”. In the Configure Reservoir menu enter a unique name, the manufacturer, the part number of the reservoir and a description, if applicable. The description can also be used for lab personal to identify the reservoir.

Dimension	Description
Reservoir Width	The outside Width (1) of the reservoir including the bottom outside flanges, i.e. the displayed dimension from back to front.
Reservoir Length	The outside Length (2) of the reservoir including the bottom outside flanges, i.e. the displayed sideward dimension.
Compartment Width	The width of the compartments, i. e. the inner sideward dimensions (3).
Compartment Length	The length of the compartments, i. e. the inner distance from back to front (4)
Compartment Count	The number of the compartments (5), e. g. 4.
Left edge offset of First Compartment	The distance between the left outside edge of the reservoir and the left edge of the first compartment, including the bottom outside flange (6).
Top edge offset of First Compartment	The distance between the top outside edge of the reservoir and the top edge of the first compartment including the bottom outside flange (7).
Distance between Compartment Positions	The distance between the compartments, from the middle of one compartment to the middle of the next (8), e. g. 9 mm. This is the compartment width plus the spacing wall.
Type	Block. Rectangular flat bottom reservoirs can be defined.
Reservoir Height	The outside height (9) of the reservoir.
Compartment Depth	The inner depth (10) of the reservoir.
max. Volume	The maximal filling volume of one compartment.

## 7.2 Plates

In the plate menu all predefined plates are listed.

**Plates**

Type: ☐ 96 wells ☐ 384 wells ☐ Other

Bottom: ☐ U-bottom ☐ F-bottom ☐ V-bottom

Volume: ☐ <1000 µl ☐ >5000 µl ☐ 1000-3000 µl ☐ 3001-5000 µl

Manufacturer: All

Part Number:

Plate	Name	Volume	Manufacturer	Part Number	Description	Favorites
	48 well F-bottom	1700 µl	GREINER	677102, 677180		★
	48 Deepwell V-Bottom Plate	7000 µl	AXYGEN CORNING	P-5ML-48-C		★
	96 well PCR Low Profile Semi Skirt	100 µl	CORNING	PCR-96-LP-AB		★
	96 Well PCR Low Profile Semi Skirt Plate	100 µl	CORNING	PCR-96-LP-AB		★
	96 Well Fast Reaction Plate	100 µl	APPLIED BIOSYSTEMS	4346907		★

To add new plates to the library, click on “Define new item”.

In the Configure Plate menu enter a unique name (e.g. well number and shape), the manufacturer and the part number of the plate.

**Configure Plate:**

Open Help

Add Plate Name:  Cancel

Save Plate

Library

Manufacturer:  Description:

Part Number:

Row Count:  8

Column Count:  12

Well Shape - Top View:  Round

Well Bottom:  V-bottom

Well Inner Diameter:  8.00 mm

Center to Center Spacing - Rows:  9.00 mm

Center to Center Spacing - Columns:  9.00 mm

Offset Row 1:  11.24 mm

Offset Column 1:  14.38 mm

Plate Height:  14.35 mm

Well Depth:  10.00 mm

Plate Width:  85.48 mm

Plate Length:  127.76 mm


max. Well Volume:  0 µl





Well V-shape Depth:  0.00 mm

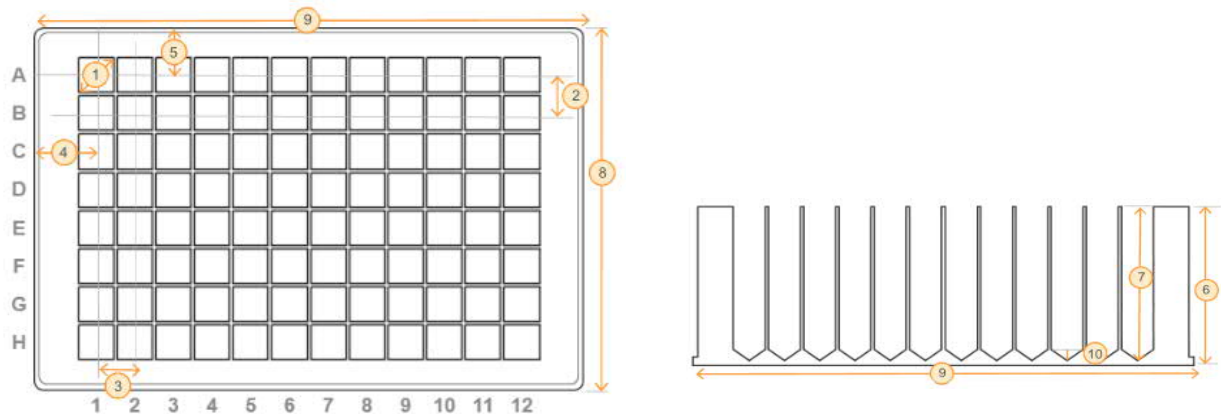

**Note:** For plates in adapters, e. g. in PCR 96 or 384 well cooling plates, define the dimensions of the inserted plate in combination with the adapter as one plate.

Enter all plate dimensions. You will be supported by live drawing of the defined plate.



Dimension	Description	
Row Count	A 96 well microplate is arranged as 8 rows by 12 columns.	
Column Count	A 384 well microplate is arranged as 16 rows by 24 columns.	
Well Shape - Top View	Shape on the top of the well, is usually round.	
Well Bottom		U-bottom = round
		F-bottom = flat
		C-bottom = curved edges of flat bottom, e. g. for easy wash
		V-bottom = conical
Well Inner Diameter	Well diameter (1) on the top of the well	
Center to Center Spacing - Rows	The distance between the rows (2) /columns (3) from the middle of one well to the middle of the next well.	
Center to Center Spacing - Columns	96 well: ANSI/SLAS-standard <sup>1</sup> : 9 mm 384 well: ANSI/SLAS-standard: 4.5 mm	
Offset Row 1	The distance between the top outside edge of the plate and the center of the first row of wells (4).	
	96 well	ANSI/SLAS-standard: 11.24 mm
	384 well	ANSI/SLAS-standard: 8.99 mm
Offset Column 1	The distance between the left outside edge of the plate and the center of the first column of wells (5).	
	96 well	ANSI/SLAS-standard: 14.38 mm
	384 well	ANSI/SLAS-standard: 12.13 mm
Plate Height	Plate height (6). ANSI/SLAS-standard: 14.35 mm ± 0.76 mm	
Well Depth	Depth in the center of the well (7).	
Plate Width	The shorter side (8), ANSI/SLAS-standard: 85.48 mm ± 0.5 mm	
Plate Length	The longer side (9), ANSI/SLAS-standard: 127.76 mm ± 0.5 mm	
max. Well Volume	The total, calculated well volume, not the working volume	
Well V-shape Depth	Depth from start of conical narrowing to well bottom only for V-bottom wells (10).	

<sup>1</sup> ANSI: American National Standards Institute, SLAS: Society for Laboratory Automation and Screening.



## 7.3 Tube Rack

To create new tube racks, select Tube Racks and click on “Define new item”. In the Configure Tube Rack menu enter a unique name, the manufacturer, the part number and a description, if applicable.

INTEGRA

Configure Tube Rack:

Open Help

Add Tube Rack Name

Cancel

Save Tube Rack

Library


96 wells





Manufacturer

Part Number

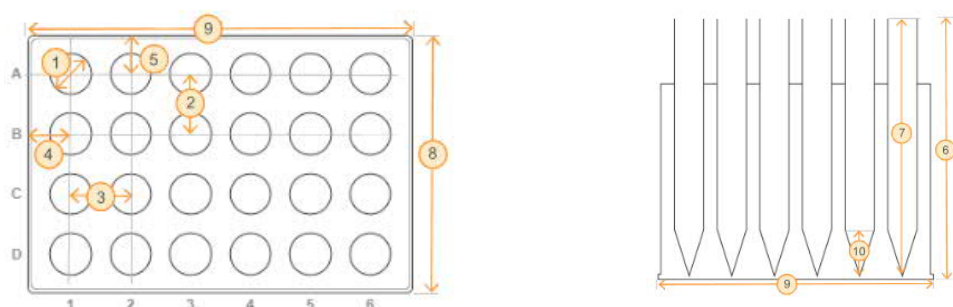Description

Enter all tube rack dimensions and the tube type. You will be supported by live drawing of the defined rack.



Dimension	Description	
Row Count	Number of rows.	
Column Count	Number of columns.	
Tube Shape - Topview	Shape on the top of the tubes, is usually round.	
Tube Shape - Bottom		U-bottom = round
		F-bottom = flat
		C-bottom = curved edges of flat bottom
		V-bottom = conical, including typical conical tube form with round bottom
Tube Inner Diameter	Tube inner diameter at the broadest position of the tube (1).	
Center to Center Spacing - Rows	The distance between the rows, from the middle of one tube to the middle of the next tube (2).	
Center to Center Spacing - Columns	The distance between the columns, from the middle of one tube to the middle of the next tube (3).	
Offset Column 1	The distance between the left outside edge of the rack and the center of the first column of tubes (4).	
Offset Row 1	The distance between the top outside edge of the rack and the center of the first row of tubes (5).	
Tube Rack Height	Total height of tube rack and tubes (6).	
Tube Depth	Depth in the center of the tube (7).	
Tube Rack Width	Width of the rack(8). ANSI/SLAS-standard: 85.48 mm ± 0.5 mm.	
Tube Rack Length	Length of the rack (9). ANSI/SLAS-standard: 127.76 mm ± 0.5 mm.	
Max. Tube Volume	The total, calculated tube volume, not the working volume.	
Tube V-shape Depth	Depth from start of conical narrowing to tube bottom, e. g. 3.1, only for V-bottom wells (10).	




## 7.4 COMBI Base

A COMBI Base is a platform where a module, i.e. a reservoir, plate or tube rack, can be placed on. It only acts as an adapter and you cannot choose a COMBI Base alone as a pipetting location (only a COMBI System see chapter 7.5)

To create a new COMBI Base, select COMBI Base and click on “Define new item”. In the Configure COMBI Base menu enter a unique name, the manufacturer, the part number of the COMBI Base and a description, if applicable.

Enter all COMBI base dimensions. You will be supported by live drawing of the defined COMBI base.

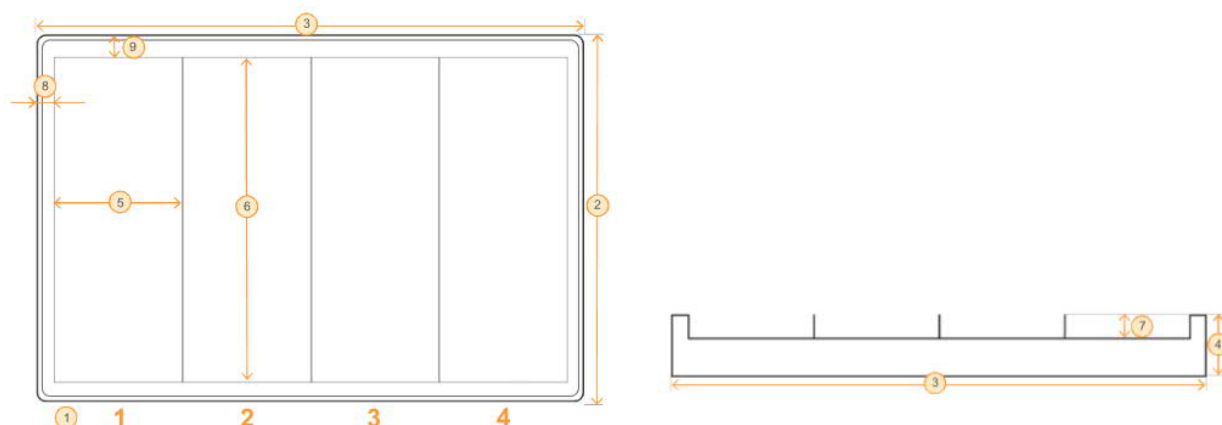

Dimension	Description
Division	<b>Regular:</b> a repeated number of identical sections are lined up. <b>Irregular:</b> sections are positioned freely within the outer dimensions of the base.
Sections	Number of sections (1).
Base Width	Width of the base (2). ANSI/SLAS-standard: 85.48 mm ± 0.5 mm.
Base Length	Length of the base (3). ANSI/SLAS-standard: 127.76 mm ± 0.5 mm.
Base Height	Total height of the base (4).
Sections Width	The width of a section (5).
Section Length	The length of a section (6).
Section Height	The height of a section (7). Negative values indicate that the bottom of the section is deeper than the whole COMBI base.
Left edge offset of section	The distance between the left outside edge of the base and the left edge of the first section, including the bottom outside flange (8).
Top edge offset of section	The distance between the top outside edge of the base and the top edge of the first section including the bottom outside flange (9).


Offset Row 1 The distance between the top outside edge of the rack and the center of the first row of tubes (5).

## 7.5 COMBI System





A COMBI System consists of a COMBI Base and individual modules, i.e. a reservoir, plate or tube rack, which are placed on the COMBI Base.

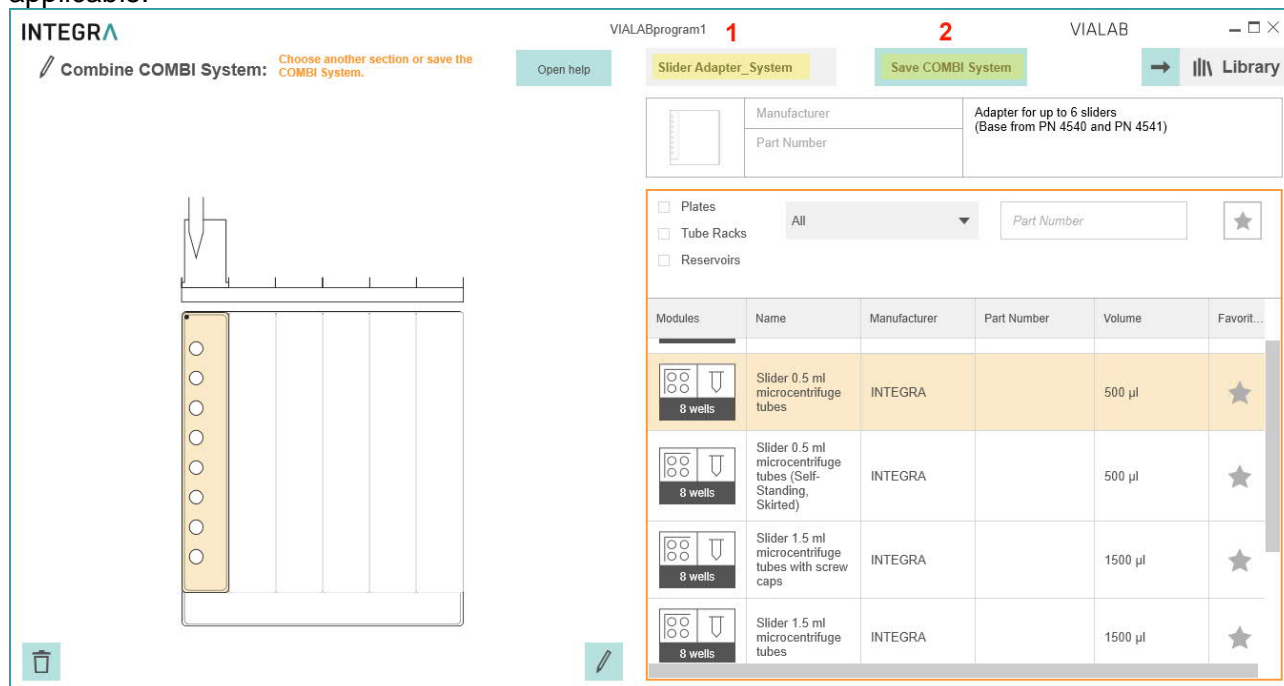To create a new COMBI System, click on “Define new item”. In the Combine COMBI System menu select a COMBI Base.

Enter a unique name, the manufacturer, the part number and a description of the COMBI System, if applicable.



The screenshot shows the INTEGRAL VIALAB software interface. At the top, there are tabs for "Slider Adapter\_System" and "Save COMBI System". Below the tabs, there is a form for entering details about the COMBI system. The form includes fields for "Manufacturer", "Part Number", and "Volume". Below these fields is a table of modules. The table has columns for "Modules", "Name", "Manufacturer", "Part Number", "Volume", and "Favorite". The table lists four modules: "Slider 0.5 ml microcentrifuge tubes", "Slider 0.5 ml microcentrifuge tubes (Self-Standing, Skirted)", "Slider 1.5 ml microcentrifuge tubes with screw caps", and "Slider 1.5 ml microcentrifuge tubes". Each module has a corresponding icon and a volume value. The "Slider 0.5 ml microcentrifuge tubes" module has a volume of 500 µl. The "Slider 0.5 ml microcentrifuge tubes (Self-Standing, Skirted)" module has a volume of 500 µl. The "Slider 1.5 ml microcentrifuge tubes with screw caps" module has a volume of 1500 µl. The "Slider 1.5 ml microcentrifuge tubes" module has a volume of 1500 µl. The interface also includes a "Combine COMBI System" button and a "Library" button.

Modules	Name	Manufacturer	Part Number	Volume	Favorit...
 8 wells	Slider 0.5 ml microcentrifuge tubes	INTEGRA		500 µl	★
 8 wells	Slider 0.5 ml microcentrifuge tubes (Self-Standing, Skirted)	INTEGRA		500 µl	★
 8 wells	Slider 1.5 ml microcentrifuge tubes with screw caps	INTEGRA		1500 µl	★
 8 wells	Slider 1.5 ml microcentrifuge tubes	INTEGRA		1500 µl	★

**1** **Name:** unique name under which the COMBI system is stored.

**2** **Save COMBI system:** save button.

Click in a section and select a Module. Continue with the next section until all sections are filled. You will be supported by live drawing of the defined COMBI System. Click on Save COMBI System.

### 7.5.1 Slanted Plate Holder

The Slanted Plate Holder (PN 4510) can accommodate plates in three different angles: 10°, 20° and 30°.

To assemble a Slanted Plate Holder with a plate select in “Material” deck position **B** or **C** on the 3 Position Universal deck. Choose the menu option “COMBI” on the right menu bar and select the COMBI Base “Slanted Plate Holder” at the desired angle.


Select either the landscape or portrait deck section on the Slanted Plate Holder and choose a compatible plate from the list on the right side.

Choose “Apply COMBI System” and use it for your program.

If you want to reuse a COMBI System, you can also create a new COMBI System in the labware library and save it with an appropriate name. Then the already assembled COMBI System (adapter and plate) can be used in your method.



**Note:** A change of the angle needs a change of the COMBI system used in the Material section. All steps in the Method section will be deleted then. We advise to define the appropriate angle in a test program first, before programming a lot of steps, which might be lost.



Combine COMBI System: Select a module.


VIALAB


→ Library

Open help

Slanted Plate Holder 20°

Apply COMBI System

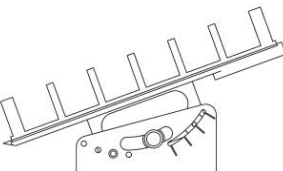



Manufacturer

4510

Slanted Plate Holder tilted at 20°

- microplates in ANSI/SLAS format can be accommodated in landscape or portrait orientation

☐ Plates






☐ Tube Racks

☐ Reservoirs

All

Part Number

★

Modules	Name	Manufacturer	Part Number	Description	Favorit...
 8 wells	8 PCR tubes				★
 24 wells	24 well F-bottom	CORNING	3337, 3524, 3526, 3527, 3473		★
 24 wells	24 well F-bottom	GREINER	662102, 662160		★
 24 wells	24 well F-bottom	TPP	92024, 92424, 92124		★
 24 wells	Back for Schott				

 **Combine COMBI System:** Select a module.

VIALABprogram1

VIALAB

Library

Manufacturer

4510

- Slanted Plate Holder tilted at 20°
- microplates in ANSI/SLAS format can be accommodated in landscape or portrait orientation

Plates







☐ Tube Racks

Reservoirs

All

Part Number



Modules	Name	Manufacturer	Part Number	Description	Favorites...
 8 wells	8 PCr tubes				
 24 wells	24 well F-bottom	CORNING	3337, 3524, 3526, 3527, 3473		
 24 wells	24 well F-bottom	GREINER	662102, 662160		
 24 wells	24 well F-bottom	TPP	92024, 92424, 92124		
 24 wells	24 well F-bottom				

*In the labware library, the Slanted Plate Holder (20°) is assembled e.g. with a 24 well plate.*

## 8 Troubleshooting

Topic	Issue	Remedy
Labware Library	I cannot place my labware item on the deck, although it is defined correctly in the labware library.	<ul style="list-style-type: none"> <li>Make sure the selected labware item is compatible with your pipette.</li> </ul>
		<ul style="list-style-type: none"> <li>Make sure your labware item fits to the selected deck position e.g. a tube rack cannot be placed on deck position A.</li> </ul>
Labware Library / Method	I created my labware item and placed it on the deck, but I cannot fill it with volume.	<ul style="list-style-type: none"> <li>Make sure the Max. Well Volume is &gt; 0 µl.</li> </ul>
Connection between pipette and computer	I cannot get a connection between pipette and VIALAB.	<ul style="list-style-type: none"> <li>Make sure the pipette is in communication mode (see chapter 6.6.1.)</li> <li>Make sure the Charging/communication stand or device is connected to the computer.</li> <li>ComModule: The Communication module for pipette and PC need to be paired if connected for the first time. Refer to Windows Help to learn how to pair Bluetooth devices. The passcode 12345 may be required.</li> <li>Make sure the drivers for the Charging/communication stand are installed. While the stand is connected, go to your System Settings -&gt; Device Manager. Check under Ports (COM &amp; LPT) if there is an entry "USB Serial Port (COM X)". If there is not, install the drivers manually.</li> </ul>
Synchronization of programs from pipette to VIALAB.	The "Synchronize" button is not activated.	<ul style="list-style-type: none"> <li>Make sure that no changes were made in VIALAB after saving the program to the pipette.</li> </ul>
Pipette status in section Transfer.	I connected the pipette and its status is "Pipette is incompatible".	<ul style="list-style-type: none"> <li>Make sure that the correct pipette model matching the current program is connected. You can only transfer programs to the pipette model used in the program.</li> </ul>

**Алматы** (7273)495-231  
**Ангарск** (3955)60-70-56  
**Архангельск** (8182)63-90-72  
**Астрахань** (8512)99-46-04  
**Барнаул** (3852)73-04-60  
**Белгород** (4722)40-23-64  
**Благовещенск** (4162)22-76-07  
**Брянск** (4832)59-03-52  
**Владивосток** (423)249-28-31  
**Владикавказ** (8672)28-90-48  
**Владимир** (4922)49-43-18  
**Волгоград** (844)278-03-48  
**Вологда** (8172)26-41-59  
**Воронеж** (473)204-51-73  
**Екатеринбург** (343)384-55-89

**Ижевск** (3412)26-03-58  
**Иваново** (4932)77-34-06  
**Иркутск** (395)279-98-46  
**Казань** (843)206-01-48  
**Калининград** (4012)72-03-81  
**Калуга** (4842)92-23-67  
**Кемерово** (3842)65-04-62  
**Киров** (8332)68-02-04  
**Коломна** (4966)23-41-49  
**Кострома** (4942)77-07-48  
**Краснодар** (861)203-40-90  
**Красноярск** (391)204-63-61  
**Курск** (4712)77-13-04  
**Курган** (3522)50-90-47  
**Липецк** (4742)52-20-81

**Киргизия** (996)312-96-26-47

**Магнитогорск** (3519)55-03-13  
**Москва** (495)268-04-70  
**Мурманск** (8152)59-64-93  
**Набережные Челны** (8552)20-53-41  
**Нижний Новгород** (831)429-08-12  
**Новокузнецк** (3843)20-46-81  
**Ноябрьск** (3496)41-32-12  
**Новосибирск** (383)227-86-73  
**Омск** (3812)21-46-40  
**Орел** (4862)44-53-42  
**Оренбург** (3532)37-68-04  
**Пенза** (8412)22-31-16  
**Петрозаводск** (8142)55-98-37  
**Псков** (8112)59-10-37  
**Пермь** (342)205-81-47

**Россия** (495)268-04-70

**Ростов-на-Дону** (863)308-18-15  
**Рязань** (4912)46-61-64  
**Самара** (846)206-03-16  
**Саранск** (8342)22-96-24  
**Санкт-Петербург** (812)309-46-40  
**Саратов** (845)249-38-78  
**Севастополь** (8692)22-31-93  
**Симферополь** (3652)67-13-56  
**Смоленск** (4812)29-41-54  
**Сочи** (862)225-72-31  
**Ставрополь** (8652)20-65-13  
**Сыктывкар** (8212)25-95-17  
**Сургут** (3462)77-98-35  
**Тамбов** (4752)50-40-97  
**Тверь** (4822)63-31-35

**Казахстан** (772)734-952-31

**Тольятти** (8482)63-91-07  
**Томск** (3822)98-41-53  
**Тула** (4872)33-79-87  
**Тюмень** (3452)66-21-18  
**Улан-Удэ** (3012)59-97-51  
**Ульяновск** (8422)24-23-59  
**Уфа** (347)229-48-12  
**Хабаровск** (4212)92-98-04  
**Чебоксары** (8352)28-53-07  
**Челябинск** (351)202-03-61  
**Череповец** (8202)49-02-64  
**Чита** (3022)38-34-83  
**Якутск** (4112)23-90-97  
**Ярославль** (4852)69-52-93